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STB-8114-P, STB-8114-G, STB-8114-C

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SONIC TOOTHBRUSH

STB-8114-P, STB-8114-G, STB-8114-C

(This manual applies to multiple color variants of the product.)

1 Intended Purpose

The product is intended only for personal oral hygiene. It is not intended for commercial use (e.g. by medical or veterinary establishments). Please follow these instructions to ensure safe use of the device.

2 Technical data

Model	STB-8114-P, STB-8114-G, STB-8114-C
Lot number	POJM180114
Power supply	1 x AAA 1.5 V --- battery
IP class	IPX6

3 Safety Instructions

1. Please read the whole manual before using the product.
2. The device is intended for household use only.
3. Use the product only for its intended purpose.
4. Do not leave the device unattended during operation.
5. Do not use any electric devices while in the bathtub or under a shower.
6. Protect the product against the penetration of water and other liquids.
7. Do not use the device if it is damaged and do not attempt to repair it.
8. The device is not a toy.
9. Keep the product out of children's reach.
10. The device should not be used by children because of the risk of electric shock.


11. Do not store the device in high humidity.
12. Do not store the product in high temperatures or sunlight.
13. Protect the device against shocks and vibration.
14. Never use strong cleaning agents, solvents or similar substances for cleaning the device.
15. Do not use the product on animals.
16. Rinse the head with water after each use.
17. Do not exert excessive pressure on your teeth or gums while brushing.
18. Do not replace the head with wet hands.
19. Remove the battery before prolonged storage.
20. Do not use any caustic agents for cleaning the product.
21. Keep this manual and, if possible, the package.
22. If the product is passed on to another user, it must be accompanied by this manual.
23. An electronic copy of this manual is available at: <http://instrukcje.vershold.com>

4 Explanation of Symbols

	Please read this manual.
	The product complies with the applicable directives of the European Union.
	Instructions for the disposal of electric and electronic devices: see the "Disposal" section.
	Product approved for contact with food
IPX6	Ingress protection class of the electrical device. The marking consists of the letters IP and two characters, of which the first denotes solid particle protection, and the second denotes liquid ingress protection. X – protection against ingress of solids not determined; 6 – full protection against water stream incoming from any direction with up to 100 L/min. flow rate.

	Trademark indicating that the manufacturer has made a financial contribution to the development and operation of a packaging materials recovery and recycling system.
	Recycling codes identifying the material from which the packaging is made – paper

5 Overview

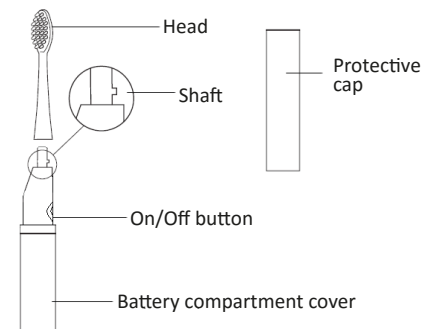6 Package Contents

Qty	Item
1	Sonic toothbrush
1	AAA 1.5 V --- battery
2	Head
1	User's Guide

Open the box and carefully remove the device. Make sure that the product kit is complete and none of its components is damaged. Check that the plastic parts are intact. If any part is damaged or missing, contact the dealer and do not use the product. Keep the package or recycle it in accordance with the applicable regulations.


Caution: Risk of suffocation! For safety, keep the packaging parts (plastic bags, cardboard boxes, Styrofoam etc.) out of children's reach.

7 Operation

7.1 Installing / replacing the battery

1. Ensure that the device is off.
2. Detach and remove the battery compartment cover.
3. Remove the old battery (if present).
4. Insert a new AAA 1.5 V --- battery, observing the polarity marked inside the compartment.

Replace the cover and tighten it, aligning the spot and the line as shown in the following illustration.

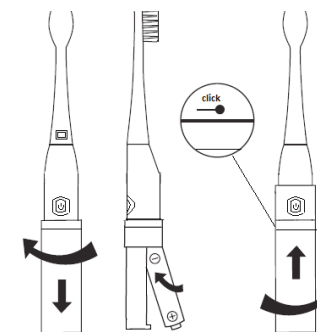Caution!

- The device is designed to be used with AAA 1,5 V --- batteries.
- The batteries are included in the set.
- When installing a new battery, observe the polarity markings in the battery compartment.
- Use the same type of batteries as recommended for this device: AAA 1.5 V --- .
- Do not dispose of used batteries with municipal waste. Put them into special containers for spent batteries.
- Do not recharge batteries not designed for that purpose (non-rechargeable batteries).
- Do not short-circuit the battery terminals.
- Never expose batteries to direct heat sources, such as bright sunlight, heaters or fire.
- Explosion hazard! Batteries must not be dismantled, put into fire or short-circuited.
- When the device is not used, always remove the batteries. This will prevent damage to the device in the event that electrolyte leaks from the batteries.
- Remove spent batteries from the device immediately.

- Spent batteries may leak, which may damage the device.
- If you touch acid from a battery, rinse your hands under running water. If any acid gets into your eyes, contact a doctor. The acid contained in the battery may cause irritation or burns.


Swallowing a battery may have fatal consequences! Keep batteries away from children and pets. Seek medical assistance immediately if a battery is swallowed.

7.2 Installing / replacing the head

- Ensure that the device is off.
- Detach the head by rotating it counter-clockwise.
- Insert and tighten a new head by rotating it clockwise.

Note:

You should feel the “click” at the end of tightening. Use only original heads to avoid damaging the toothbrush.

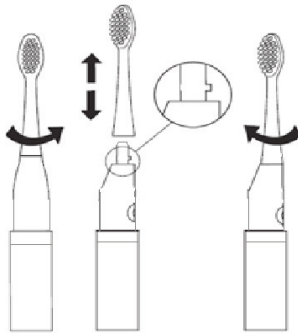7.3 Brushing

- Wet the toothbrush head.
- Apply toothpaste to the brush.
- To avoid splashing toothpaste around, insert the brush in your mouth before turning it on.
- Turn the device on by pressing the On/Off button.
- Start brushing from the gums on the extreme left or right of your dental arches.
- Brush all teeth one after another, starting from the gum.
- Switch to the inner side and repeat these steps.
- Then scrub the teeth chewing surfaces carefully.
- Once you are finished brushing, turn off the device with the On/Off button.
- Rinse the head under running water and dry it.

Notes:

The recommended time of application of the electric toothbrush is 3 minutes.

At first, your gums may bleed slightly as a result of the new type of irritation. This is normal, but if this condition persists for more than one week, consult your dentist.

8 Cleaning and Maintenance

General rules

Correct and regular cleaning will improve the safety and extend the lifespan of the product.

Caution: Turn the device off and leave it to cool down before cleaning to avoid electric shock or burn.

Caution: Do not use any strong chemical, alkaline, abrasive or disinfecting agents for cleaning the device as they could damage its surface.

- Clean the toothbrush head thoroughly after brushing by rinsing it under running water.
- For optimum hygiene, you should clean the toothbrush head with a denture cleaning agent every 2-4 weeks.

9 Service

The device does not contain any user-serviceable parts. Do not attempt to make any repairs. Just call the Service.

10 Handling and storage

Clean the device according to the instructions.

Keep the device in a dry and well-ventilated place inaccessible to children.

Protect the device against vibration and shocks during handling.

11 Disposal

All the packaging materials are recyclable and labelled as such. Dispose of the packaging in accordance with the local regulations.

Keep the materials out of children’s reach, as the materials can pose a hazard.

Proper disposal:

- According to the WEEE Directive (2012/19/UE), the crossed-out wheellie bin symbol (shown on the left) is used to label all electric and electronic devices requiring segregation.
- Do not dispose of the spent product with domestic waste: hand it over to a electric and electronic device collection and recycling point. The crossed-out

wheellie bin symbol placed on the product, instruction manual or package communicates this requirement.

- Plastics contained in the device can be recycled in accordance with their specific marking. By recycling materials and spent equipment you help protect the environment.
- Information on electric and electronic devices collection locations is available from local government agencies or from the dealer.
- Spent or fully discharged single-use and rechargeable batteries must be discarded in dedicated labelled bins, handed over to hazardous waste collectors or returned to electric equipment dealers.

12 CE Declaration of Conformity

The device has been designed, manufactured and marketed in compliance with the EMC Directive (2014/30/EU). Therefore, the product bears the “CE” mark and has a declaration of conformity, which is available to market supervision agencies for inspection.

13 Warranty

In order to claim a product, you must deliver the product to the Customer Service Center at any Biedronka store.

If you have any product-related questions or problems, write us at:

infolinia@vershold.com

or call us at: **+48 667 090 903**.

Your opinion is important to us. Please rate our product on our page:
www.vershold.com/opinie

- The Manufacturer of this product grants a 24-month warranty valid from the purchase date. Should you find any defect, submit a complaint at the place of purchase.
- A defective/damaged product means a product failing to provide the features described in the Instruction Manual due to intrinsic properties of the appliance.
- The Beneficiary under the warranty shall deliver the appliance to the Customer Service Point at a relevant store in order to submit a complaint about the equipment. When asserting your rights under the warranty, you need to present proof of purchase (purchase receipt, invoice) together with a description of the device defect in as much detail as possible, in particular of the external manifestations of the defect.

4. The warranty excludes products with physical damage other than that caused by the manufacturer or distributor, including specifically any damage resulting from incorrect usage or force majeure.

5. The warranty does not cover products from which the tamper-proof seal has been removed and products to which repairs, alterations or structural modifications have been attempted.

6. The warranty does not cover components subject to natural wear and tear during operation.

7. It is recommended that the complete product be returned in order to facilitate defect verification.

8. This warranty for the sold consumer product does not exclude, restrict or suspend the rights of the buyer related to the regulations concerning implied warranty for defects pursuant to the Civil Code of 23 April 1964 (Journal of Laws Dz.U.2014.121).

9. The warranty shall be applicable within the territory of the Republic of Poland.

Manufacturer (“Warrantor”):

VERSHOLD POLAND Sp. z o.o.

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02-092 Warszawa, Poland

Made in China