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
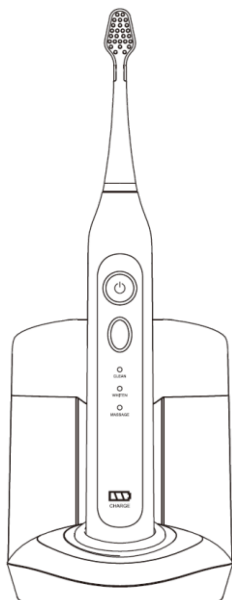
User manual

Rechargeable Sonic Power Toothbrush with UV Sanitizer

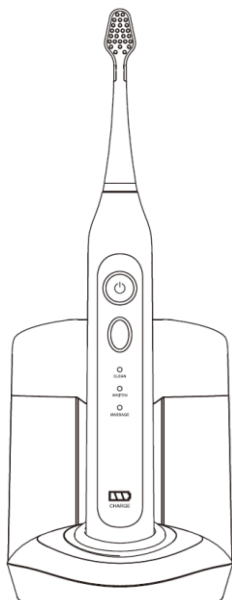Specification and colours may slightly differ from the product shown in the illustrations.

Table of Contents

1. INTRODUCTION.....	16
2. INTENDED USE.....	16
3. TECHNICAL DATA.....	16
4. PRINCIPLES OF SAFE USE.....	17
5. EXPLANATIONS OF THE SYMBOLS.....	19
6. CONTENTS OF THE SET.....	20
7. CONSTRUCTION.....	21
8. USE.....	22
9. CLEANING AND MAINTENANCE.....	26
10. PROBLEMS.....	27
11. STORAGE AND TRANSPORT.....	28
12. WARRANTY AND SERVICE.....	28
13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS.....	28
14. MANUFACTURER.....	28

1. INTRODUCTION

Thank you for choosing our product. We are confident that its use will be a great pleasure for you. Please read the user manual carefully before using the product, in particular the instructions on safe use. Please keep this manual in order to be able to refer to it while using the product.

2. INTENDED USE

The sonic toothbrush is intended exclusively for personal oral cavity care in the household. Do not use the sonic toothbrush commercially, industrially, or in medical or veterinary practice.

3. TECHNICAL DATA

Rechargeable Sonic Power Toothbrush with UV Sanitizer

Model	AL-5202
Lot number	202008
Number of pulsations	40 000 micro vibrations per minute
UV-C wave length	254 nm
Number of programs	33
Program duration	2 minutes with 30 seconds interval
Use time after full charge	ca. 10 days (10 days x 2 minutes x 2 cleaning modes a day)
Ingress protection (toothbrush)	IPX7
Material	ABS
Charging base power supply	100-240 V~ 50/60 Hz, 3 W, power cord 1.3 m
Battery	Ni-MH 3,6 V \approx , 700 mAh
Charging type	wireless
Charge time	~24 h
Accessories	3 interchangeable tips, charging base
Dimensions of charging base	98 x 100 x 127 mm
Dimensions of toothbrush (with the tip):	248 x 32 mm
Weight of toothbrush (with the tip):	140g

4. PRINCIPLES OF SAFE USE

Read the entire user manual carefully before use.

WARNING! DANGER OF INJURY! RISK GROUP 3

The product emits UV radiation! Avoid exposing your eyes and skin to UV light with the product uncovered.

This appliance may be used by children of at least 8 years of age and people with reduced physical and mental capabilities, or lack of experience and knowledge of the appliance, if supervised; or instructed about safe use so that the hazards involved are understood.

Unattended children must not clean or maintain the unit.

Children must not play with this device.

- The device emits UV-C ultraviolet radiation with a wavelength of 254 nm.
- Avoid exposure of the eyes and skin to ultraviolet (UV) radiation.
- Use this device only when it is properly assembled (also applies to the charging base).
- Protect people and animals! UV radiation may cause injury to eyes and skin! Never risk exposing yourself, other people or animals to direct or reflected UV radiation. In case of exposure to radiation, contact your doctor immediately.
- Use the device only as described in this user manual!
- Disconnect the power plug from the outlet before installing, cleaning, or maintaining the unit.
- Do not use a damaged device. Do not repair yourself if damaged. Take the device to a professional for repair or deliver it to the place of purchase.
- Do not use this device in locations where it is prohibited to use electrical equipment, or where it may cause harmful interference or other hazards.
- Do not switch on or use the device close to flammable materials.
- Do not use the machine to clean surfaces other than teeth, gums and

tongue.

- Keep the device out of the reach of children.
- The manufacturer is not responsible for any problems caused by unauthorized modifications to the equipment.
- When using the device, ensure that there are no objects left near the place where the toothbrush will be placed that could be damaged, broken or knocked over. The power supply base cable should be routed so that it does not come into contact with objects that generate high temperatures. It should also not be subjected to crushing, pulling, or other damage. The charging base cable should not interfere with the movement of persons and pets, as this can lead to tripping, tipping over or injury.
- Avoid exposing the device to very high/low temperatures (below 0°C/32°F and over 40°C/104°F). Plastic parts may deform under extreme temperatures, which may affect the operation of the device or lead to complete damage.
- Do not expose the charging base and its power cord to water or other liquids. This may adversely affect the operation and lifetime of equipment components or cause a short circuit or even electric shock. Do not immerse the toothbrush in water.
- Avoid exposing the battery to very high/low temperatures (below 0°C/32°F and over 40°C/113°F). Extreme temperatures may influence capacity and battery life. Avoid contact with liquids and metal objects as it may cause partial or complete battery damage. Use the battery only for its intended purpose. Do not destroy, impair or throw the battery into flames - it might be dangerous and may cause fire. Worn-out or impaired battery should be placed in a dedicated container. Never open the battery. Please dispose of battery according to the instructions. Disconnect an unused charging base from power grid.

You can obtain a free printed version of the manual with the large font size or download it from the web address below:

<https://mptech.eu/qr/hoffen>

5. EXPLANATIONS OF THE SYMBOLS

	<p>The device is marked with a crossed-out rubbish bin, in accordance with the European Directive 2012/19/WE on used electrical and electronic Equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled. Collecting this type of waste in designated locations, and the actual process of their recovery contribute to the protection of natural resources. Proper disposal of used electrical and electronic equipment is beneficial to human health and environment. To obtain information on where and how to dispose of used electrical and electronic equipment in an environmentally friendly manner user should contact their local government office, the waste collection point, or point of sale, where the equipment was purchased.</p>			
	<p>In accordance with directive 2013/56/UE on the disposal of batteries and accumulators, this product is marked with the symbol of a crossed-out waste bin. The symbol indicates that the batteries used in this product should not be disposed with regular household waste, but treated according to the legal directives and local regulations. It is not allowed to dispose of batteries and accumulators with unsorted municipal waste. Users of batteries and accumulators must use the available collection points of these items, which allows them to return, recycle and disposal. Within the EU, the collection and recycling of batteries and accumulators is a subject to separate procedures. To learn more about existing procedures of recycling of batteries and accumulators please contact your local office or an institution for waste disposal or landfill.</p>			
	<p>The mark placed on the product indicates that the product complies with the basic requirements and other provisions of the so-called "Directive of the 'New Approach' Directive of the European Union.</p>			
				<p>The packaging of the device can be recycled in accordance with local regulations.</p>

	Symbol for cardboard (packaging material).
	For indoor use.
	Symbol for Class II equipment. Equipment where protection against electric shock is provided not only by basic insulation but also by additional safety measures such as double insulation or reinforced insulation.
<p data-bbox="86 619 191 667">IPX7</p>	The standard for immersion in water for up to 30 minutes at a depth of up to 1 meter (applies to the lowest part of the device). Items that meet this standard can be used in shallow water swimming, bathing and showering. The device may seep when immersed in water more than 1 metre deep, when immersed for more than 30 minutes, or when subjected to a stream of water under very high pressure.
<p data-bbox="86 783 191 831">IPX4</p>	Degree of protection against water: is resistant to water splashes from any direction.

6. CONTENTS OF THE SET

After unpacking, check that all listed parts are in the package and that the device is working properly. Included in the set:

1. Sonic toothbrush
2. Interchangeable tips (3 pcs)
3. Charging base with power cord
4. USER MANUAL

If any of the above items is missing or damaged, please contact your dealer.



We advise to keep the container in case of any complaints. Keep the packaging away from children - risk of swallowing and/or suffocation.

Replaceable brush available for purchase from the manufacturer.

Contact: (+48 71) 71 77 400 or pomoc@myphone.pl

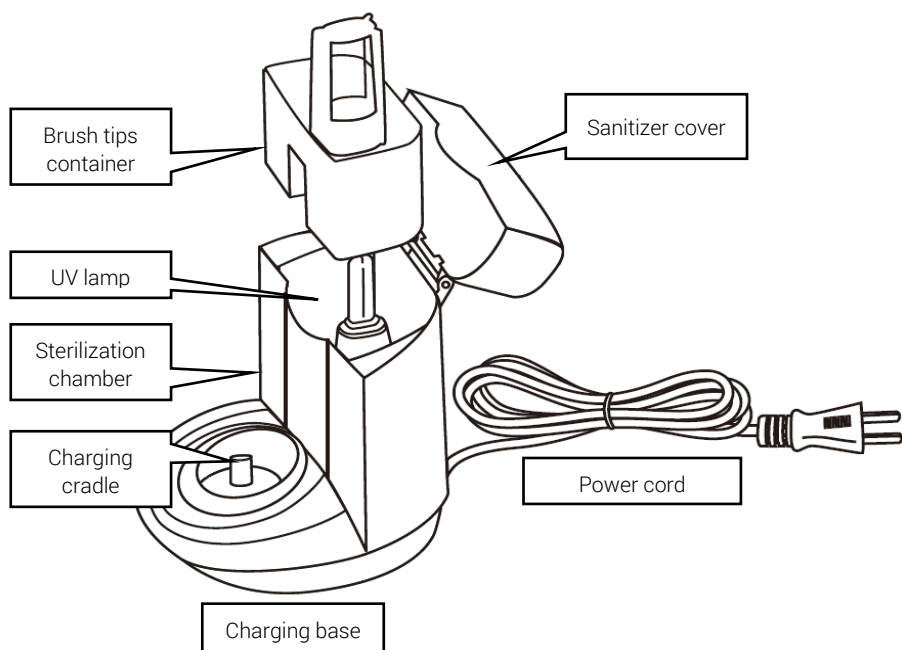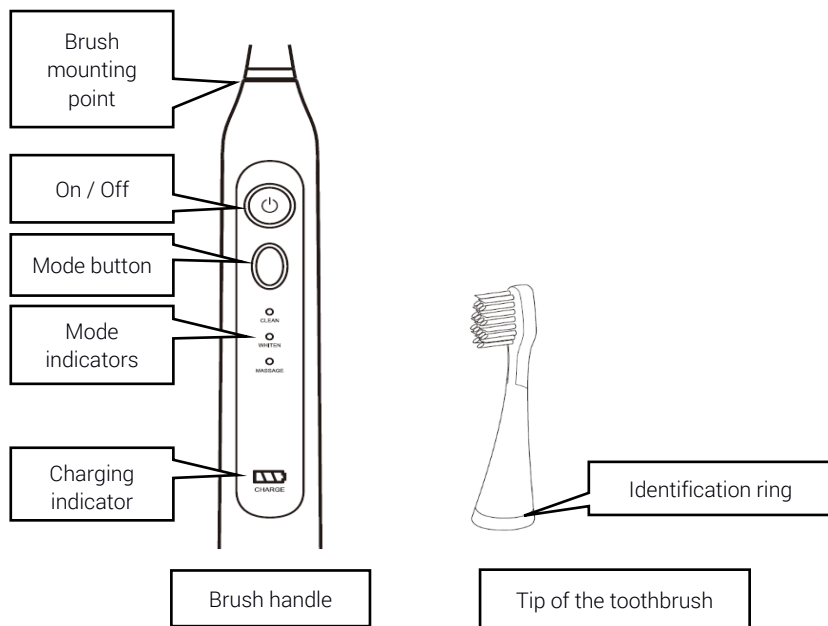7. CONSTRUCTION

THE GENERAL APPEARANCE OF THE DEVICE


8. USE

-  *Mounting and dismounting of the kit components is only permitted when the unit is switched off.*
-  *Carry out all the steps gently, careful not to damage the kit components.*
-  *Before and after usage, check each component for signs of damage. Damaged components must not be used, as this may cause damage to the user's body.*
-  *Do not use the tip for more than 3 months after first use. Teeth or gums might get injured due to excessive use of single brush or bristles wearing off.*

MOUNTING AND DISMOUNTING OF BRUSH TIPS

Take out the toothbrush and one of the tips from the packaging. Insert the tip onto the protruding spindle from the brush handle so that the bristles of the tips point in the same direction as the switch and diodes on the brush handle. Press the two parts together in a firm motion. This process is illustrated in the following figure. A small space of approx. 0.5 mm will be left between the brush tip and the brush handle. This is the result of the design of the device. Do not be press against the brush handle forcing the tip to come into contact with the brush.

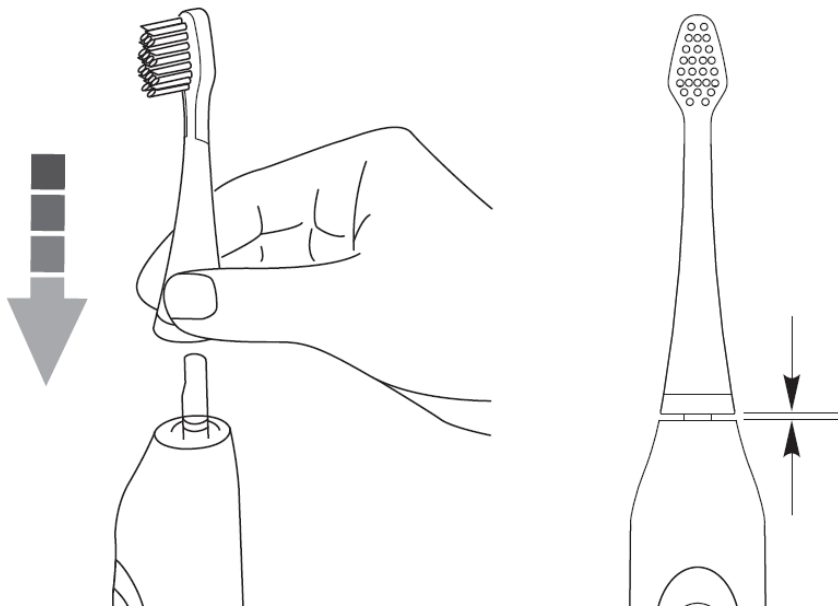
To remove the tip, gently pull the tip upwards until it is detached from the brush handle.



CHARGING THE TOOTHBRUSH BATTERY

Keep your hands dry before operating the charging base and charging the battery in the device.

Before charging the battery for the first time, remove the charging base from the package. Unwind the charging base cable gently and place it on a stable horizontal surface so that the cable reaches the power socket. Plug the base cable plug into a power outlet.

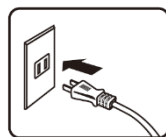Place the bottom of the brush handle in the charging cradle. Depending on the current battery charge level, the charge light can take one of the four possible states shown below:

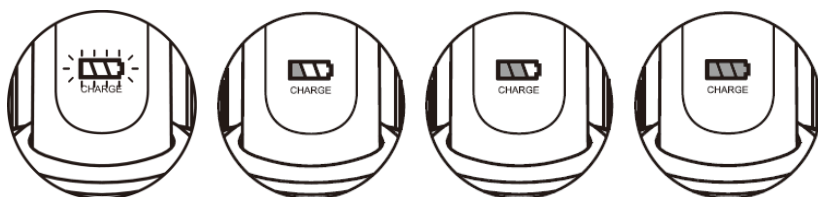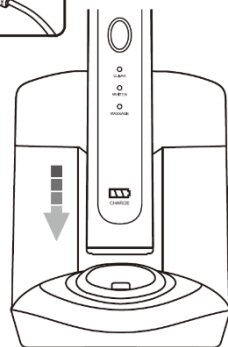- Three blinking segments of the charging light indicate a charge of less than 1/3 of the battery capacity.
- One constantly illuminated segment and two blinking segments indicate a charge between 1/3 and 2/3 of the battery capacity.
- Two constantly illuminated segments and one flashing indicates a charge above 2/3 of the battery capacity.
- Three permanently illuminated segments of the charging light indicate that the battery is fully charged.

Low energy level in the battery is indicated by flashing of the leftmost segment of charging indicator. If the battery is completely discharged, the device may not turn on or stop immediately after turning it on.

Do not charge or store the unit below 0°C or above 40°C, otherwise the quality and length of the charge will be affected.

The device may become slightly warm during charging; which is normal.

Try not to discharge the battery completely, as this can lead to permanent damage to the battery and brush handle.

TOOTHBRUSH MODES

The toothbrush offers three modes of operation:

- **Clean** – Standard mode for daily brushing.
- **White** – mode facilitating teeth whitening and discoloration removal, brightening, front teeth polishing. Recommended for people who drink discolouring drinks and smoke cigarettes.
- **Massage** – intensive gum massage after a standard brushing session, designed for healthy gums.


The toothbrush has a built-in timer and after every 30 seconds it informs you with a short break that the time has passed. The standard washing session provides 4 such intervals, which gives a total of 2 minutes of washing/massage.


TECHNIQUE OF CLEANING WITH A SONIC BRUSH

To use the sonic toothbrush effectively, follow these instructions.


Apply a small dose of your favourite toothpaste to the bristles of the brush. In the case of sonic toothbrushes, the amount of toothpaste may be smaller than in the case of regular toothbrushes.

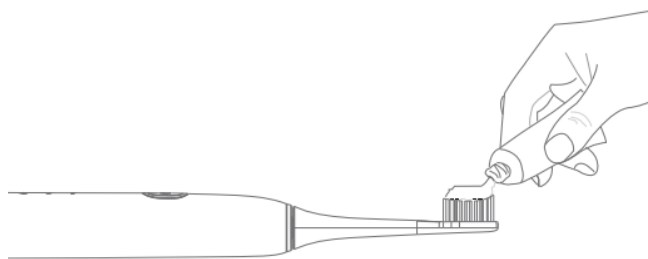The bristles of the toothbrush should be directed at an angle of 45 degrees so that the bristles have the ability to clean the side surfaces of the tooth and gums.

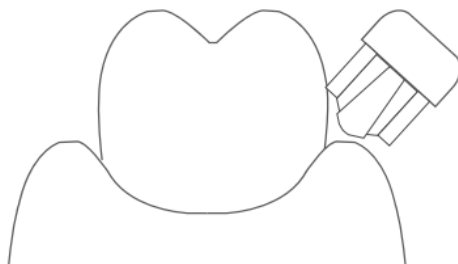Switch on the device and guide the toothbrush tip slowly, pressing it gently against the tooth surface and touching gums. Single cleaning cycle lasts 4x30s. It is worth dividing the surface of your teeth into four parts so that you can clean each one for 30 seconds. After 2 minutes of operation, the brush will turn off.

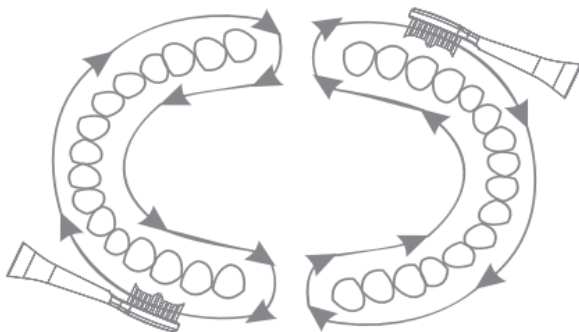
 *During the first dozen or so uses of the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. In such case, it is recommended to use "Massage" mode. If, after 2 weeks of using the toothbrush in the recommended way gums still bleed, stop using the device and consult a dentist.*

USING THE STERILIZATION CHAMBER

 *Sterilize only clean tips.*

The charging base is equipped with a sterilization chamber which enables sterilization of brush tips with UV-C radiation after each use.

Before using this function, gently shake the tip to get rid of any excess water remaining after cleaning. Open the sterilization chamber cover, put the tips in the tip container, placing them on the pins at the bottom of the container so that the brushes face the UV lamp. Make sure the charging base is plugged to power supply. Close the sterilization chamber to start the sanitizing process. At this point a glow with a slightly blue tint will be visible from under the cover of the sterilization chamber, signalling the operation of the UV lamp. The sterilization process takes about 7 minutes, after that time the sterilizer will automatically switch off and the glow will no longer be visible.

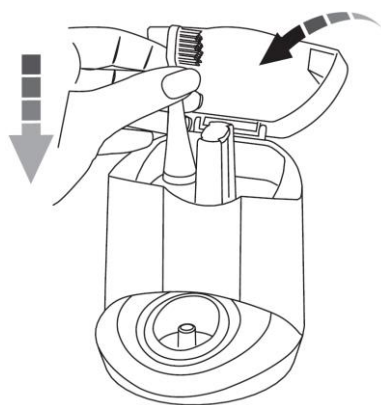In well-lit rooms, glow can be difficult to see.

Sterilize only the tips supplied with the device or compatible, bought from the manufacturer. Do not place other objects in the sterilizing chamber.

The sterilizer will automatically switch on every time its flap is closed.

9. CLEANING AND MAINTENANCE

CLEANING OF THE EQUIPMENT COMPONENTS

Regular and proper cleaning and maintenance will ensure a long service life and trouble-free operation of the equipment.

Cleaning and maintenance may only be carried out when the appliance is switched off.


It is not allowed to clean any part of the toothbrush in the dishwasher.

Do not immerse the brush handle under water.

The unit's components can be cleaned under running water at low pressure and temperature. In the case of heavy soiling, mild detergents, such as non-concentrated dishwashing liquids, can be used. Then thoroughly dry the surfaces of the components with a paper or textile towel and allow to dry completely in an airy place. The brushes should be washed under running water and then allowed to dry. From time to time you can wash the tip using a mouthwash.

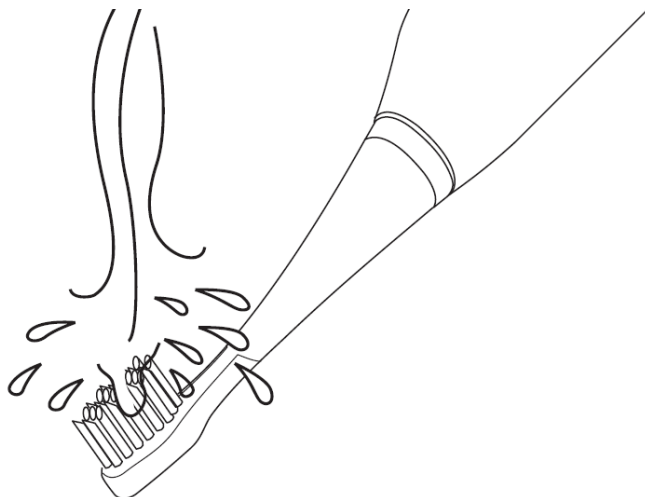10. PROBLEMS

In case of any problems with the device, please refer to the following tips.

The toothbrush does not switch on	Check if the battery is charged.
Vibrations are too weak	Check that you are not pressing the brush against the teeth too hard. Check whether the brush tip or brush housing are damaged. Check that the tip is not accidentally too firmly seated on the brush's pin, which may hinder its movement.
Despite the toothbrush has been placed in a charging base, the brush does not charge.	Check that the cable or plug coming out of the charging base is not damaged. Check if the charging base is correctly connected to the power supply. Perhaps the service life of the battery is over - contact the manufacturer's service department.
The sanitizer does not switch on	Check if the charging base is correctly connected to the power supply. Check that the sterilization chamber cover is closed. If the UV lamp is damaged, please contact the manufacturer's service department.
The gums bleed after using the toothbrush.	Check that you do not press the tip too hard against your gums. Check that you do not press the tip against your gums at a wrong angle. <i>During the first dozen or so uses of the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. In this case, it is recommended to use "Massage" mode. If, after 2 weeks of using the toothbrush in the recommended way gums still bleed, stop using the device and consult a specialist dentist.</i>
If the problem persists, contact the manufacturer's service department.	

11. STORAGE AND TRANSPORT

STORAGE OF EQUIPMENT

If you are not going to use the device for a long time, after disconnecting the charging base from the power source, cleaning and drying its components, store it in a dry room away from heat and moisture sources. It is best to place the device in its original packaging, away from children and pets.

TRANSPORT OF THE DEVICE

It is recommended to transport the device in its original packaging. Do not expose the unit to excessive vibration or shock during transportation. During operation, only move the unit when it has been disconnected from the power supply. Do not expose the device to high temperatures or moisture.

12. WARRANTY AND SERVICE

The product (the toothbrush and the charging base) is covered by a 3 year warranty from the date of sale of the product. In case of a defect, the product should be returned to the shop where it was purchased, remembering to provide the proof of purchase. The product should be complete and in its original packaging. Consumables (replaceable brushes) are not covered by the manufacturer's warranty.

13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS

 mPTech Sp. z o.o. hereby declares that the product complies with all relevant provisions of the so called the 'New Approach' Directive of the European Union.

14. MANUFACTURER

mPTech Sp. z o.o.
Nowogrodzka street 31
00-511 Warsaw
Poland

Helpline: (+48 71) 71 77 400
Manufactured in China
Lot number: **202008**