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The device is available in various colour versions. The colours of the product may vary slightly from the illustrations in this manual.

TABLE OF CONTENTS

1. INTRODUCTION11
 2. INTENDED USE11
 3. TECHNICAL DATA11
 4. PRINCIPLES OF SAFE USE11
 5. EXPLANATIONS OF THE SYMBOLS13
 6. THE CONTENTS OF THE SET14
 7. CONSTRUCTION14
 8. USE15
 9. CLEANING AND MAINTENANCE17
 10. PROBLEMS18
 11. STORAGE AND TRANSPORT19
 12. WARRANTY AND SERVICE19
 13. DECLARATION OF CONFORMITY WITH THE NEW APPROACH DIRECTIVES19
 14. PRODUCENT19

1. INTRODUCTION

Thank you for choosing our product. We are confident that its use will be a great pleasure for you. Please read the user manual carefully before using the product, in particular the instructions on safe use. Please keep this manual in order to be able to refer to it while using the product.

2. INTENDED USE

The toothbrush is intended exclusively for personal teeth care in the household. Do not use toothbrush commercially, industrially, or in medical or veterinary practice.

3. TECHNICAL DATA

Device: Oscillating toothbrush for children
 Model: AL-1392
 Lot number: 202005
 Number of pulsations: 800/min
 Tightness degree: IPX5 (protection from water)
 Tips: 3 interchangeable
 Power supply: 3V (2 x AA battery) 1.5

4. PRINCIPLES OF SAFE USE

Read these guidelines carefully. Not following them may be dangerous or illegal.

IMPORTANT SAFETY INSTRUCTIONS
 The device should not be used and maintained by children under 8 years of age and persons with limited mental or physical predispositions, unless under supervision or after handing over the instructions for use of the device in a safe manner, understanding the danger posed by using the device. This device is not a toy. Cleaning and maintenance must not be carried out by unsupervised children.

DO NOT RISK
 The device is intended for domestic use and is not suitable for commercial, industrial, medical, or veterinary purposes. Do not use this device in locations where it is prohibited to use electrical equipment, or where it may cause harmful interference or other hazards. Do not switch on or use the device close to flammable materials. Do not use the machine to clean surfaces other than teeth, gums and tongue. Keep the device out of the reach of children.

SURROUNDINGS

When using the device, ensure that there are no objects left near the place where the toothbrush will be placed that could be damaged, broken or knocked over.

QUALIFIED SERVICE

Only qualified personnel or an authorized service center may install or repair this product. Repairing by an unauthorized or unqualified service may cause device damage and warranty loss. The manufacturer is not responsible for problems caused by unauthorized modifications to the equipment made by the user.

BATTERY AND ACCESSORIES

Avoid exposing the battery to very high/low temperatures (below 0°C/32°F and over 40°C/104°F). Extreme temperatures can affect service life of the batteries. Avoid contact with liquids and metal objects as it may cause partial or complete battery damage.

Use the battery only according to its intended purpose. Do not destroy, impair or throw the battery into flames - it might be dangerous and may cause fire. Place worn or damaged batteries in a special container. Never open the battery. Dispose of the battery according to the instructions. Use the original manufacturer's accessories only.

5. EXPLANATIONS OF THE SYMBOLS

The device is marked with a crossed-out garbage bin, in accordance with the European Directive 2012/19/EU on used electrical and electronic equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled. Collecting this type of waste in designated locations, and the actual process of their recovery contribute to the protection of natural resources. Proper disposal of used electrical and electronic equipment is beneficial to human health and environment. To obtain information on where and how to dispose of used electrical and electronic equipment in an environmentally friendly manner user should contact their local government office, the waste collection point, or point of sale, where the equipment was purchased.

In accordance with directive 2013/56/EU on the disposal of batteries and accumulators, this product is marked with the symbol of a crossed-out waste bin. The symbol indicates that the batteries used in this product should not be disposed with regular household waste, but treated according to the legal directives and local regulations. It is not allowed to dispose of batteries and accumulators with unsorted municipal waste. Users of batteries and accumulators must use the available collection points of these items, which allows them to return, recycle and disposal. Within the EU, the collection and recycling of batteries and accumulators is a subject to separate procedures. To learn more about existing procedures of recycling of batteries and accumulators please contact your local office or an institution for waste disposal or landfill.

The mark placed on the product indicates that the product complies with the basic requirements and other provisions of the so-called Directive of the European Union.

The packaging of the device can be recycled in accordance with local regulations.

Symbol for cardboard (packaging material).

The OTHER symbol indicates plastics other than PET, HDPE, PVC, LDPE, PP and PS.

Degree of protection against a stream of water (water poured onto the housing from any direction does not cause harmful effects)

6. THE CONTENTS OF THE SET


After unpacking, check that all replaced parts are in the package and that the device is working properly. Included in the set:

- 1. Oscillating Toothbrush for Children
- 2. 3 x replaceable brushes
- 3. 2 x AA 1,5 V battery
- 4. User manual

If any of the above items is missing or damaged, please contact your dealer. We advise to keep the container in case of any complaints. Keep the packaging away from children - risk of swallowing and/or suffocation.

7. CONSTRUCTION

Device appearance

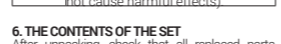8. USE

 ⚠️ Mounting and dismounting of the kit components is only permitted when the unit is switched off.
 ⚠️ Carry out all the steps gently, careful not to damage the kit components.
 ⚠️ Before and after work, check each component for signs of damage. Damaged components must not be used, as this may cause damage to the user's body.

CHANGING BRUSH TIPS

⚠️ Do not use the tip for more than 3 months from the first use. Teeth or gums might get injured due to excessive use of single brush or bristles wearing off. The toothbrush in the package has one tip already mounted. If it wears out, replace it. To replace the tip, rotate the tip against the handle by 90° and then separate the tip from the handle by pulling it in opposite directions. This process is illustrated in the illustration on the right.

To insert a new tip, rotate it 90° from the handle and point its hole at the stem protruding from the toothbrush holder. Insert it to the end and twist the hair forward.

INSERTING / REPLACING BATTERIES

⚠️ There are 2 AA batteries included in the kit.

⚠️ Make sure your hands are dry before inserting and replacing batteries in the device.

⚠️ Improper insertion of batteries will only result in the device not working properly.

CAUTION!

- Do not charge non-rechargeable batteries (other than accumulators).
- Do not short-circuit the battery power terminals.
- Use batteries of the same type as recommended for this device.
- When installing new batteries, remember to check

the +/- polarity of the batteries.
 • Remove discharged batteries immediately from the unit, as they may leak and cause damage.
 • Always remove the batteries when not in use - this will prevent the product from being damaged if the electrolyte leaks from the battery.
 • Always replace the entire battery pack with new ones.
 • Keep batteries away from children and animals.
 • Do not dispose of your used batteries in a waste bin, but in special containers for used batteries and accumulators.
 • Never use used batteries in combination with new ones.
 • Avoid exposing the battery to very high/low temperatures (below 0°C/32°F and over 40°C/104°F). Extreme temperatures can affect their service life.
 • Avoid contact with liquids and metal objects, as this may result in total or partial damage to the battery.
 • Never expose batteries to direct sunlight, radiators, fire, excessive heat etc.
 • Never dispose of batteries into a fire - they may explode.
 • Never immerse batteries in water. Do not touch batteries with wet hands.
 • If hands come into contact with acid from the battery, rinse hands under running water. If electrolyte gets into your eyes, contact your doctor. Acid contained in batteries may cause skin irritation or chemical burns.
 • Swallowing batteries can be fatal! Keep batteries away from children and pets. Immediately seek medical attention if the battery is swallowed.
CAUTION! In order to ensure adequate safety measures and operating conditions, the battery should only be installed by adults. Before the first use, a new set of two 1.5 V AA batteries must be inserted into the battery compartment.
 To remove the battery compartment cover, hold the battery compartment cover firmly with fingers of one hand, and the brush holder with the other hand.
 When replacing the batteries, remove the complete set of used batteries and replace them with 2 new 1.5V AA batteries . After inserting or replacing the battery, insert the battery compartment cover

into the brush holder until it clicks into place.

TEETH BRUSHING

⚠️ The toothbrush is not able to clean any areas between teeth, so it is recommended to use floss after brushing - once a day.
CAUTION! During the first dozen or so uses of the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. If, after 2 weeks of using the toothbrush in the recommended way gums still bleed, stop using the device and consult a specialist dentist.

⚠️ If a child uses the brush, make sure it spits out all the toothpaste leftovers after washing.

To switch the unit on and off, use the buttons shown in one of the previous illustrations. To use your toothbrush effectively, follow these instructions. Apply a small dose of your favourite toothpaste to the bristles of the brush. In the case of oscillating toothbrushes, the amount of toothpaste may be smaller than in the case of ordinary toothbrushes. A portion the size of a pea is usually sufficient.

The bristles of the toothbrush should be directed at an angle of 45 degrees so that the bristles have the ability to clean the side surfaces of the tooth, gums and interdental spaces. Guide the toothbrush tip slowly, pressing it gently against the tooth surface and joining it with the gums, and make gentle sweeping movements from the gums to the crown of the tooth. Then clean the contact surfaces of the upper and lower jaw teeth. The usual cleaning time should be about 2 minutes.

9. CLEANING AND MAINTENANCE

CLEANING OF THE EQUIPMENT COMPONENTS

⚠️ Regular and proper cleaning and maintenance will ensure a long service life and trouble-free operation of the equipment.

⚠️ Cleaning and maintenance may only be carried out when the appliance is switched off.

⚠️ It is not allowed to clean any part of the

toothbrush in the dishwasher.
 The unit can be cleaned under running water at low pressure and temperature. In the case of heavy soiling, mild detergents, such as non-concentrated dishwashing liquids, can be used. Then dry the surface of the device thoroughly with a paper or textile towel and leave it to dry completely in a well-ventilated place. The brushes should be washed under running water and then allowed to dry. From time to time you can wash the tip using a mouthwash.

10. PROBLEMS

In case of any problems with the device, please refer to the following tips.
The toothbrush does not switch on
 Check the batteries. If they depleted, replace them as described in this manual.
 If the battery is replaced by a new one, the toothbrush does not work, check for mechanical damage or flooding.

Vibrations are too delicate

Check that you are not pressing the brush against the teeth too hard.
 Check whether the brush tip or brush housing are damaged.
 Make sure that the tip is not too firmly seated on the pin in the brush (the distance between the tip and brush is too small for the tip to move freely).

The gums bleed after using the toothbrush

Check that you do not press the tip too hard against your gums.
 During approximately first dozen of uses the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. If, after 2 weeks of using the toothbrush in the recommended way gums still bleed, stop using the device and consult a specialist dentist.

If the problem persists, contact the manufacturer's service department.

11. STORAGE AND TRANSPORT
STORAGE OF EQUIPMENT
 If you are not going to use the device for a long time, after removing the battery, cleaning and drying its parts, store it in a dry room away from heat and moisture sources. It is best to place the device in its original packaging, away from children and pets.

TRANSPORT OF THE DEVICE

It is recommended to transport the device in its original packaging. Do not expose the unit to excessive vibration or shock during transportation. When handling, only carry the unit after it has been switched off, holding it securely by the housing. Do not expose the device to heat or moisture.

12. WARRANTY AND SERVICE

The product (the sonic toothbrush) is covered by a 3 year warranty from the date of sale of the product. In case of a defect, the product should be returned to the shop where it was purchased, remembering to provide the proof of purchase. The product should be complete and preferably delivered in its original packaging. Consumables (replaceable brushes and batteries) are not covered by the manufacturer's warranty.

13. DECLARATION OF CONFORMITY WITH THE NEW APPROACH DIRECTIVES

mPTech Ltd. hereby declares that the product complies with all relevant provisions of the so-called the 'New Approach' Directive of the European Union.

14. MANUFACTURER

mPTech Sp. z o.o. Helpline: (+48 71) 71 77 400
 Nowogrodzka 31 street Manufactured in China
 00-511 Warsaw Lot number: 202005
 Poland

You can obtain a free printed version of the manual with the large font size or download it from the web address below:
www.mptech.eu/gr/hoffen/choose-product/

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