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USER MANUAL RECHARGABLE OSCILLATING TOOTHBRUSH

| Table of Contents | |
|--|----|
| 1. INTRODUCTION..... | 14 |
| 2. INTENDED USE..... | 14 |
| 3. TECHNICAL DATA..... | 14 |
| 4. PRINCIPLES OF SAFE USE..... | 15 |
| 5. EXPLANATIONS OF THE SYMBOLS..... | 16 |
| 6. CONTENTS OF THE SET..... | 18 |
| 7. CONSTRUCTION..... | 18 |
| 8. USE..... | 19 |
| 9. CLEANING AND MAINTENANCE..... | 21 |
| 10. TROUBLESHOOTING..... | 22 |
| 11. STORAGE AND TRANSPORT..... | 22 |
| 12. WARRANTY AND SERVICE..... | 23 |
| 13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS..... | 23 |
| 14. MANUFACTURER..... | 23 |

1. INTRODUCTION

Thank you for choosing our product. We are confident that its use will be a pleasure for you. Please read the user manual carefully before using the product, in particular the instructions on safe use. Please keep this manual in order to be able to refer to it while using the product.

2. INTENDED USE

The toothbrush is intended exclusively for personal oral cavity care in the household. Do not use the product commercially, industrially, or in medical or veterinary practice.

3. TECHNICAL DATA

| Rechargeable oscillating toothbrush | |
|-------------------------------------|---|
| Model | AL-6202 |
| Lot number | 202011 |
| Velocity | 8800 oscillations/min. |
| Number of programs | 3 |
| Program duration | 2 minutes with 30 seconds interval |
| Use time after full charge | ca. 10 days (10 days x 2 minutes x 2 cleaning modes a day) |
| Ingress protection | IPX7 |
| Material | ABS |
| Cradle charger power supply | USB 5 V = 200 mA |
| Built-in battery | Ni-MH 2.4 V = 600 mAh |
| Charging type | wireless |
| Charge time | ~24 h |
| Accessories | 2 interchangeable tips, charging base |
| Weight of toothbrush with the tip: | 110g |

4. PRINCIPLES OF SAFE USE

Read the entire user manual carefully before use.

This appliance may be used by children of at least 8 years of age and people with reduced physical and mental capabilities, or lack of experience and knowledge of the appliance, if supervised; or instructed about safe use so that the hazards involved are understood.

Unattended children must not clean or maintain the unit.

Children must not play with this device.

The device is intended for use only with the included charging base. Do not use other charging bases to charge the device.

- Use this device only when it is properly assembled (also applies to the charging base).
- Use the device only as described in this user manual!
- Disconnect the power plug from the outlet before cleaning, or maintaining the unit.
- Do not use a damaged device. Do not repair by yourself if damaged. Take the device to a professional for repair or deliver it to the place of purchase.
- Do not use this device in locations where it is prohibited to use electrical equipment, or where it may cause harmful interference or other hazards.
- Do not use the machine to clean surfaces other than teeth, gums and tongue.
- Keep the device out of the reach of children.
- The manufacturer is not responsible for any damages caused by unauthorized modifications to the equipment.
- When using the device, ensure that there are no objects left near the place where the toothbrush will be placed that could be damaged,

broken or knocked over. The power supply base cable should be routed so that it does not come into contact with objects that generate high temperatures. It should also not be subjected to crushing, pulling, or other damage. The charging base cable should not interfere with the movement of persons and pets, as this can lead to tripping, tipping over or injury.

- Avoid exposing the device to very high/low temperatures (below 0°C/32°F and over 40°C/104°F). Plastic parts may deform under extreme temperatures, which may affect the operation of the device or lead to complete damage.
- Do not expose the charging base and its power cord to water or other liquids. This may adversely affect the operation and lifetime of equipment components or cause a short circuit or even electric shock. Do not immerse the toothbrush in water.
- Avoid exposing the battery to very high/low temperatures (below 0°C/32°F and over 40°C/104°F). Extreme temperatures may influence capacity and battery life. Avoid contact with liquids and metal objects as it may cause battery damage. Use the battery only for its intended purpose. Do not destroy, impair or throw the battery into flames - it might be dangerous and may cause fire. Worn-out or impaired battery should be placed in a dedicated container. Never open the battery. Please dispose of battery according to the instructions. Disconnect an unused cradle charger from power grid.

You can obtain a free printed version of the manual with the large font size or download it from the web address below:

<https://mpotech.eu/qr/hoffen>

5. EXPLANATIONS OF THE SYMBOLS

The device is marked with a crossed-out rubbish bin, in accordance with the European Directive 2012/19/WE on used electrical and electronic Equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled. Collecting this type of waste in designated locations, and the actual process of their recovery contribute to the protection of natural resources. Proper disposal of used electrical and electronic equipment is

| | | |
|---|--|--|
| | <p>beneficial to human health and environment. To obtain information on where and how to dispose of used electrical and electronic equipment in an environmentally friendly manner user should contact their local government office, the waste collection point, or point of sale, where the equipment was purchased.</p> | |
|  | <p>In accordance with directive 2013/56/UE on the disposal of batteries and accumulators, this product is marked with the symbol of a crossed-out waste bin. The symbol indicates that the batteries used in this product should not be disposed with regular household waste, but treated according to the legal directives and local regulations. It is not allowed to dispose of batteries and accumulators with unsorted municipal waste. Users of batteries and accumulators must use the available collection points of these items, which allows them to return, recycle and disposal. Within the EU, the collection and recycling of batteries and accumulators is a subject to separate procedures. To learn more about existing procedures of recycling of batteries and accumulators please contact your local office or an institution for waste disposal or landfill.</p> | |
|  | <p>The mark placed on the product indicates that the product complies with the basic requirements and other provisions of the so-called "Directive of the 'New Approach' Directive of the European Union.</p> | |
|  |  |  <p>The packaging of the device can be recycled in accordance with local regulations.</p> |
|  | <p>Symbol for cardboard (packaging material).</p> | |
|  | <p>Symbol for PVC - Polyvinyl chloride (packaging material).</p> | |
| <p>IPX7</p> | <p>The standard for immersion in water for up to 30 minutes at a depth of up to 1 meter (applies to the lowest part of the device). Items that meet this standard can be used in shallow water swimming, bathing and showering. The device may seep when immersed in water more than 1 metre deep, when immersed for more than 30 minutes, or when subjected to a stream of water under very high pressure.</p> | |

6. CONTENTS OF THE SET

After unpacking, check that all listed parts are in the package and that the device is working properly. Included in the set:

1. Rechargeable oscillating toothbrush
2. Interchangeable tips (2 pcs)
3. Charging base with USB power cord
4. User manual

If any of the above items is missing or damaged, please contact your dealer.

We advise to keep the container in case of any complaints. Keep the packaging away from children - risk of swallowing and/or suffocation.

Replaceable brush available for purchase from the manufacturer.

Contact: (+48 71) 71 77 400 or pomoc@myphone.pl

7. CONSTRUCTION

THE GENERAL APPEARANCE OF THE DEVICE

8. USE

Mounting and dismounting of the kit components is only permitted when the unit is switched off.

Carry out all the steps gently, careful not to damage the kit components.

Before and after work, check each component for signs of damage. Damaged components must not be used, as this may cause damage to the user's body.

Do not use the tip for more than 3 months after first use. Teeth or gums might get injured due to excessive use of single brush or bristles wearing off.

MOUNTING AND DISMOUNTING OF BRUSH TIPS

Take out the toothbrush and one of the tips from the packaging. Turn the tip 45° in the vertical axis from the front of the brush holder. Place the tip in the mounting location and then turn 45° to the left.

To remove the tip, turn it 45° clockwise and then gently pull it upwards until it is detached from the brush holder.

CHARGING THE TOOTHBRUSH BATTERY

Keep your hands dry before operating the charging base and charging the battery in the device.

Before charging the battery for the first time, remove the charging base from the package. Unwind the charging base cable gently and place it on a stable horizontal surface so that the cable reaches the power socket. Plug the base cable plug into a power outlet (using USB adapter).

Place main part of the tooth brush in the charging base. Depending on the current battery charge level, the charge light can take one of the two possible states: blue or red.. Blue light means that the battery is sufficiently charged for operation, red light means that the battery needs to be recharged. The battery in the brush holder is fully charged when the charging light goes out.

The device indicates a low energy level in the battery even when it is switched on - the charging light will then light up red. In rare cases, the charging light will light up blue as soon as it is switched on, but will change to red after a while. In the event of extreme battery discharge, the device may not turn on or might turn off in the middle of a cleaning programme, not completing it.

Do not charge or store the unit below 0°C or above 40°C, otherwise the quality and length of the charge will be affected.

The device may become slightly hot during charging; which is normal.

Try not to discharge the battery completely, as this can lead to permanent damage to the battery and brush.

The built-in battery is not user replaceable. To replace the battery, contact the manufacturer's service department.

TOOTHBRUSH MODES

The toothbrush offers three modes of operation:

- **Clean** – Standard mode for daily brushing.
- **Sensitive** – the more gentle oscillation mode, designed for people with delicate gums.
- **Massage** – intensive gum massage after a standard brushing session, designed for healthy gums.

Pressing the switch button again immediately after turning on the device will switch between modes. When all modes have been activated, pressing the switch again will switch the brush off.


The toothbrush has a built-in timer and after every 30 seconds it informs you with a short break that the time has passed. The standard washing session provides 4 such intervals, which gives a total of 2 minutes of washing/massage.

TECHNIQUE OF CLEANING WITH THE OSCILLATING BRUSH


To use your toothbrush effectively, follow these instructions.

Apply a small dose of your favourite toothpaste to the bristles of the brush. In the case of oscillating toothbrushes, the amount of toothpaste may be smaller than in the case of ordinary toothbrushes.

The bristles of the toothbrush should be directed at an angle of 45 degrees so that the bristles have the ability to clean the side surfaces of the tooth and gums.

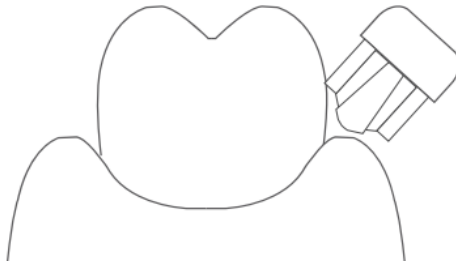Switch on the device and guide the toothbrush tip slowly, pressing it gently against the tooth surface and touching gums. Single cleaning cycle lasts 4x30s. It is worth dividing the surface of your teeth into four parts so that you can clean each one for 30 seconds. After 2 minutes of operation, the brush will turn off.

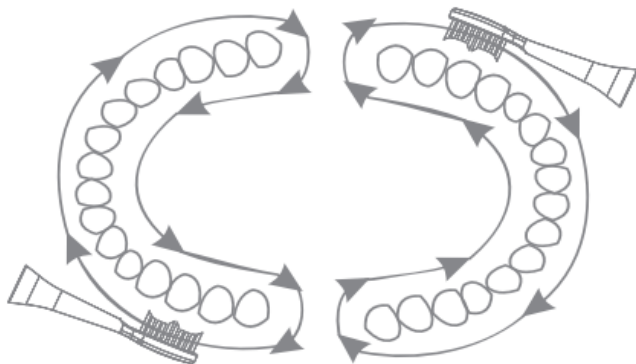 *During the first dozen or so uses of the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. In this case, it is recommended to use "Massage" mode. If, after 2 weeks of using the toothbrush, in the recommended way, gums still bleed, stop using the device and consult a specialist dentist.*

9. CLEANING AND MAINTENANCE

CLEANING OF THE EQUIPMENT COMPONENTS

 *Regular and proper cleaning and maintenance will ensure a long service life and trouble-free operation of the equipment.*

 *Cleaning and maintenance may only be carried out when the appliance is switched off.*

 *It is not allowed to clean any part of the toothbrush in the dishwasher.*

 *Do not immerse the brush holder under water.*

The unit's components can be cleaned under running water at low pressure and temperature. In the case of heavy soiling, mild detergents, such as non-concentrated dishwashing liquids, can be used. Then thoroughly dry the surfaces of the components with a paper or textile towel and allow to dry completely in an airy place. The brushes should be washed under running water and then allowed to dry. From time to time you can wash the tip using a mouthwash.

10. TROUBLESHOOTING

In case of any problems with the device, please refer to the following tips.

| | |
|--|--|
| <p>The toothbrush does not switch on</p> | <p>Check that the battery is charged.</p> |
| <p>The oscillations are too delicate</p> | <p>Check that you are not pressing the brush against the teeth too hard. Check whether the brush tip or brush housing are damaged. Check that the battery is charged properly.</p> |
| <p>Despite the toothbrush has been placed in a charging base, the brush does not charge.</p> | <p>Check that the cable or plug coming out of the charging base is not damaged. Check if the charging base is correctly connected to the power supply. Perhaps the service life of the battery is over - contact the manufacturer's service department.</p> |
| <p>The gums bleed after using the toothbrush.</p> | <p>Check that you do not press the tip too hard against your gums. Check that you do not press the tip against your gums at a bad angle. <i>During the first dozen or so uses of the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. In this case, it is recommended to use "Massage" mode. If, after 2 weeks of using the toothbrush in the recommended way gums still bleed, stop using the device and consult a specialist dentist.</i></p> |
| <p>If the problem persists, contact the manufacturer's service department.</p> | |

11. STORAGE AND TRANSPORT

STORAGE OF EQUIPMENT

If you are not going to use the device for a long time, after disconnecting the charging base from the power source, cleaning and drying its components, store it in a dry room away from heat and moisture sources. It is best to place the device in its original packaging, away from children and pets.

TRANSPORT OF THE DEVICE

It is recommended to transport the device in its original packaging. Do not expose the unit to excessive vibration or shock during transportation. During

operation, only move the unit when it has been disconnected from the power supply. Do not expose the device to high temperatures or moisture.

12. WARRANTY AND SERVICE

The product (the toothbrush and the charging base) is covered by a 3 year warranty from the date of sale of the product. In case of a defect, the product should be returned to the shop where it was purchased, remembering to provide the proof of purchase. The product should be complete and in its original packaging. Consumables (replaceable brushes) are not covered by the manufacturer's warranty.

13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS

mPTech sp. z o.o. hereby declares that the product complies with all relevant provisions of the so-called the 'New Approach' Directive of the European Union.

14. MANUFACTURER

mPTech Sp. z o.o.
Nowogrodzka 31 street
00-511 Warsaw
Poland

Helpline: (+48 71) 71 77 400
E-mail: pomoc@mpotech.eu
Manufactured in China
Lot number: **202011**