58.8mm

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85mm

HOFFEN

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Model: ST-1257-W, ST-1257-M, ST-1257-P

CE IPX7

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HOFFEN

SONIC TOOTHBRUSH
Model: ST-1257-W, ST-1257-M,
ST-1257-P

The manual applies to multiple colour variants of the product

1. INTENDED PURPOSE

The sonic toothbrush is intended only for personal oral hygiene. It can be used by adults and children over 3 years of age. Please follow these instructions to ensure safe use of the device

2. SPECIFICATIONS

Lot number	POJM210257
Model	ST-1257-W, ST-1257-M, ST-1257-P
Power supply	1 x AAA 1,5 V --- battery (included)
IP class	IPX7
Net weight:	51,4 g
Gross weight:	93,4 g

3. SAFETY INSTRUCTIONS


1. Please read the whole manual before using the product.
2. The device is intended for household use only.
3. It is recommended that children not be left unattended when using the toothbrush.
4. Use the product only for its intended purpose.
5. Do not use the device if it is damaged, and do not attempt to repair it.
6. The device is not a toy.
7. Do not store the device in high humidity.
8. Do not store the device in high temperatures or sunlight.
9. Protect the device against shocks and impacts.
10. Never use strong cleaning agents, solvents or similar substances to clean the device.
11. Do not use the product on animals.
12. Clean the toothbrush head after every use.
13. Do not exert excessive pressure on teeth or gums while brushing.
14. Do not replace the head when your hands are wet.
15. Remove the battery before prolonged storage.
16. Keep this manual and, if possible, the package.
17. If the product is passed on to another user, it must be accompanied by this manual.
18. An electronic copy of this manual is available at: <http://instrukcje.vershold.com>

4. EXPLANATION OF SYMBOLS

	Please read this manual.
	The product complies with the applicable directives of the European Union.
	Instructions for the disposal of electric and electronic devices: see the "DISPOSAL" section.
	Protection against a short-time immersing in water under defined conditions of pressure and time (30 minutes up to 1 meter of submersion).
	Trademark indicating that the manufacturer has contributed financially to the construction and operation of a packaging materials recovery and recycling system.

	Recycling code identifying the material the packaging is made of – cardboard, polyethylene terephthalate.
	Segregation sign - Labeling indicating the need for waste segregation.
	Plastic/metal waste - Designation of the container to which the waste should go.
	Paper waste - Designation of the container to which the waste should go.

5. OVERVIEW

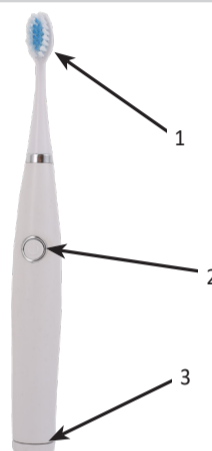1. Head
2. On/off button
3. Battery compartment cover

6. PACKAGE CONTENTS


Qty	Item
1	Sonic toothbrush
1	AAA 1,5 V --- battery
3	Head (1 installed, 2 extra)
1	User's Guide

Open the box and carefully remove the device. Make sure that the product kit is complete and none of its components is damaged. Check that the plastic parts are intact. If any part is damaged or missing, contact the dealer and do not use the product. Keep the package or recycle it in accordance with the applicable regulations.

Attention! For the sake of children's safety, do not leave any packaging parts accessible and unattended (plastic bags, cardboard boxes, Styrofoam etc.). Risk of suffocation!

7. OPERATION

- ### 7.1 Installing / replacing the battery
1. Ensure that the device is turned off.
 2. Rotate the battery compartment cover, as shown in the illustration below.
 3. Remove the battery compartment cover.
 4. Remove the old battery (if present).
 5. Insert a new AAA 1.5 V --- battery, observing the polarity marks in the lower part of the compartment.
 6. Replace the cover and tighten it, as shown in the illustration below.

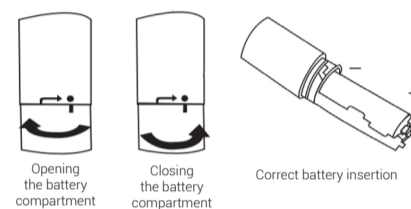- Caution!**
- The device is designed to be used with AAA 1.5 V --- batteries.
 - The batteries are included.
 - When replacing batteries, observe the polarity marked on the top of the battery compartment cover.
 - Use only the recommended batteries (1,5 V, AAA ---).
 - Do not dispose of old batteries with domestic waste. Discard them in a dedicated bin available at many stores.
 - Never re-charge single-use (non-rechargeable) batteries.
 - Do not short-circuit the battery terminals or contacts.

8. CLEANING AND MAINTENANCE


- Always remove the batteries before prolonged storage of the device to prevent damage from possible electrolyte leak.
- Remove spent batteries from the device without delay.
- Spent batteries can leak, which can damage the product.
- After exposure of hands to leaked electrolyte, immediately wash the hands under running water. In the event that battery acid comes into contact with eyes, immediately seek

medical attention. The acid can irritate or burn eyes.

- Swallowing a battery can be fatal! Keep batteries out of the reach of children and pets. If you suspect your child or pet has swallowed a battery, seek medical attention immediately.

7.2 Installing / replacing the head

1. Ensure that the device is turned off.
2. Detach the head by rotating it clockwise.
3. Install the new head.

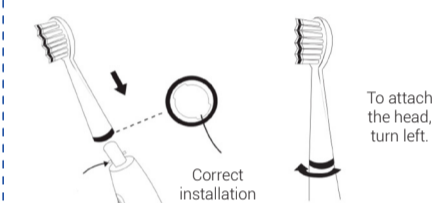58.8mm

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Caution! You should feel the "click" at the end of tightening. Use only original heads to avoid damaging the toothbrush. The set contains 3 heads in total.

7.3 Brushing

- Wet the toothbrush head.
- Apply some toothpaste to the wet brush.
- To avoid splashing toothpaste around, insert the brush in your mouth before turning it on.
- Turn the device on by pressing the On/Off button. The brush will vibrate.
- Start by brushing the outside surfaces of the teeth. Press the brush head slightly against the gums and brush the teeth starting from the gums.
- Brush all teeth one after another,

starting from the gum.

- Switch to the inner side and repeat these steps.
- Then brush the backs of your molars carefully.
- Once you are finished brushing, turn off the device with the On/Off button. The brush will turn off.
- After you have completed brushing the teeth, rinse the head under running water and dry it.

Notes:
The recommended length of use of the electric toothbrush is 3 minutes. At first, your gums may bleed slightly as a result of the new type of irritation. This is normal. However, if this condition persists for more than one week, consult your dentist.

8. CLEANING AND MAINTENANCE

- General rules**
Correct and regular cleaning will improve the safety and extend the lifespan of the product.
- Warning!** Turn the device off and leave it to cool down before cleaning to avoid electric shock or burn.
- Caution!** Do not use any strong chemical, alkaline, abrasive or disinfecting agents to clean the device as they could damage its surface.
- Clean the toothbrush head thoroughly after brushing by rinsing it under running water.
 - For the sake of hygiene, it is recommended that the brush head be changed every 3 months.

9. SERVICE

The device does not contain any user-serviceable parts. Do not attempt to make any repairs. Just call the Service.

10. HANDLING AND STORAGE

Clean the device according to the instructions. Keep the device in a dry and well-ventilated place inaccessible to children. Protect the device against vibration and shocks during transportation.

11. DISPOSAL

All the packaging materials are recyclable and labelled as such. Dispose of the packaging in accordance with local regulations. Keep the packaging materials out of children's reach, as the materials can pose a hazard.

Proper disposal:

1. According to the WEEE Directive 2012/19/EU, the crossed-out wheellie bin symbol (shown on the left) is used to label all electric and electronic devices requiring segregation.
2. Do not dispose of the spent product with domestic waste: hand it over to an electric and electronic device collection and recycling point. The crossed-out wheellie bin symbol placed on the product, instruction manual or package communicates this requirement.
3. Plastics contained in the device can be recycled in accordance with their specific marking. By recycling materials and spent equipment you help protect the

4. Information on electric and electronic devices collection locations is available from local government agencies or from the dealer.
5. Spent or fully discharged single-use and rechargeable batteries must be discarded in dedicated labelled bins, handed over to hazardous waste collectors or returned to electric equipment dealers.
6. Failure to comply with the recycling rules may lead to negative effects on the environment and human health.

12. CE DECLARATION OF CONFORMITY

The device has been designed, manufactured and marketed in compliance with the EMC Directive. Therefore, the product bears the "CE" mark and has a declaration of conformity, which is available to market supervision agencies for inspection.

13. WARRANTY

To submit a complaint about the device, deliver it to a Customer Service Point at any Biedronka store.

If you have any product-related questions or problems, write us at: infolinia@vershold.com or call us at: +48 667 090 903

Your opinion is important to us. Please rate our product on our page: www.vershold.com/opinie

1. The Warrantor of these products grants a 24-month warranty valid from the date of purchase. Should you find any defect, submit a complaint at the place of purchase.
2. A defective/damaged product means a product failing to provide the features described in the User's Manual due to intrinsic properties of the device.
3. The Beneficiary under the Warranty shall deliver the appliance to the Customer Service Point at a given store in order to submit a complaint regarding the device. When exercising your warranty rights, present proof of purchase (receipt, invoice) and describe the defect in as much

4. The warranty excludes products with physical damage other than that caused by the manufacturer or distributor, including specifically any damage resulting from incorrect usage or force majeure.
5. The warranty does not cover products whose protective seal has been removed or to which repairs, modifications or structural changes have been made or attempted.
6. The warranty does not cover components subject to natural wear and tear during operation.

7. It is recommended that the entire product be turned over upon submission of a complaint, so as to facilitate service operations.
8. This warranty for the sold consumer product does not exclude, restrict or suspend any rights of the Buyer stemming from warranty regulations covering defects in goods sold pursuant to the Civil Code of 23 April 1964 (Journal of Laws Dz.U.2014.121).
9. This Warranty is valid on the territory of Poland..

Manufacturer (Warrantor)
VERSHOLD POLAND Sp. z o.o.
ul. Żwirki i Wigury 16A
02-092 Warszawa, Poland
Made in China