# HOFFEN

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## 1. INTENDED USE

Skin scrubber painlessly and thoroughly cleanses the skin, accelerates its nourishment and regeneration.

Cavitation – a phenomenon occurring while applying ultrasound waves to wet skin. It created microscopic gas bubbles, followed by a fast pressure change and several other processes which result in breaking down any dead cells in the stratum corneum of the epidermis.

Apply the cavity peeling treatment to your face, neck and neckline area.

Follow this manual to safely install and use the device.

This device is intended for indoor use only. It may not be used for commercial purposes.

## 2. TECHNICAL SPECIFICATION

Lot number	POJM210311
Model	SS-1311
Battery parameters	3.7 V  ; 500 mAh (built-in, non-replaceable)
USB charging cable parameters	5 V  ; 500 mA
Charging time	Approx. 130 minutes

## 3. SAFETY INSTRUCTIONS

1. Read the entire manual before use. Following the instructions ensures safe installation and use of the device.
2. A device for private home use. Do not use for other purposes not in accordance with its intended purpose.
3. The device is not a toy. Make sure that children do not play with the device.
4. When using electrical appliances, especially when children are present, basic safety precautions should be observed to avoid the risk of fire.
5. Cleaning and maintenance of the device must not be performed by children.
6. Do not immerse the device in water or other liquids while cleaning or operating.
7. Never keep the device under running water. Do not use the device while bathing or showering.

8. **Caution!** Contraindications for application:
  - pregnancy and breast feeding
  - cancer
  - hypertension
  - heart failure
  - pacemaker
  - epilepsy
  - fever
  - bacterial, viral and fungal infections of the skin
  - purulent skin inflammation, dermatitis
  - Roaccutane acne treatment – 3 months break
  - metal implants in treated area
  - orthodontic braces
9. Avoid prolonged and intensive use of the device in one spot to prevent irritation and redness of treated skin.
10. If skin irritation occurs, stop the treatment immediately.
11. Do not use on damaged or irritated skin.
12. Keep the instructions and, if possible, the packaging.
13. If we pass the device on to a third party, we must pass these operating instructions along with it.
14. The electronic version of the manual is available at:  
**<http://instrukcje.vershold.com/>**
15. **Caution!** Battery Warning:
  - Do not put the used batteries in a rubbish bin, but in special containers for spent single-use and rechargeable batteries.
  - Never expose the battery to direct heat sources, such as excessive sunlight, heaters or fire.
  - Danger of explosion! Batteries must not be disassembled, thrown into fire or short-circuited.
  - Do not disassemble, open, or damage the battery housing.
  - Do not subject the battery to mechanical shocks.
  - If any acid from the battery comes into contact with your hands, rinse them under running water. If any acid gets into your eyes, contact a doctor. Battery acid can cause irritation or burns.

## 4. EXPLANATION OF SYMBOLS

Read the manual

The product complies with the applicable directives of the European Union

Disposal of waste electrical and electronic equipment – see the DISPOSAL section of this manual.

Used, spent batteries and rechargeable batteries must be discarded in specially marked containers, handed over to special waste collection points or to electrical equipment dealers

The device is intended for indoor use only.

Trademark indicating that the manufacturer has made a financial contribution to the development and operation of a packaging materials recovery and recycling system.

Designation of the packaging material – cardboard

Segregation sign – Labeling indicating the need for waste segregation.


Paper waste – Designation of the container to which the waste should go.

Designation of type of waste bin for the packaging – PLASTIC/METAL.

Symbols in monochrome printout only for illustration.

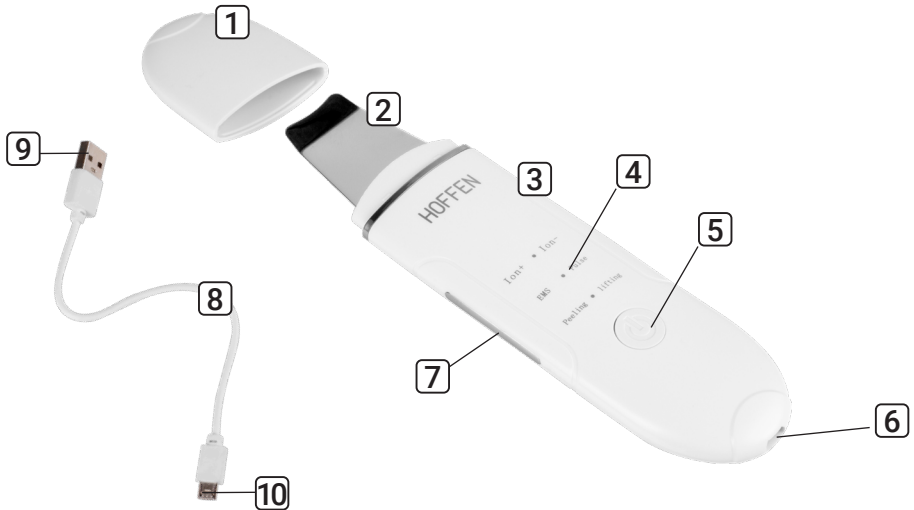1. Cover
2. Spatula
3. Device together with built-in battery
4. Function display:
  - Peeling / lifting (cleansing / lifting)
  - EMS / pulse (EMS massage)
  - Ion + (positive ionization)
  - Ion - (negative ionization)
5. On / off / function switch
6. USB charging socket
7. Ion sensor
8. USB cable
9. USB end

## 6. SET CONTENTS

The device set should contain the following items:

- 1 x Cavity peeling machine with cap
- 1 x USB charging cable
- 1 x User manual

Open the packaging and carefully take out the device. Ensure that the package is complete and none of its contents are damaged. Check that the plastic parts are not broken. If any part is damaged or missing, contact the dealer and do not use the product. Keep the package or recycle it in accordance with local regulations.

**Caution! For safety, keep the packaging parts (plastic bags, cardboard boxes, Styrofoam etc.) out of the reach of children. Risk of suffocation!**

## 7. USE

### 7.1 Charging

The device signals a low battery level with a double beep.

To charge the device, connect its charging cable to:

a) Computer

- Connect the charging cable plug to the charging port of the device, then insert the USB connector to your computer's USB port. Make sure that the computer is turned on.
- After connecting the cable to the computer, the LED will light up and all of the LEDs will flash.
- When the device is fully charged, the LEDs stop flashing and the LED indicator lights off.

b) Power adapter (not supplied)

- Connect the charging cable plug to the charging port of the device, then insert the USB connector into a USB port of a power adaptor with parameters compliant with the device label. Connect the power adapter to a mains outlet with the same power characteristics as those stated

on the power adapter label.

- After connecting the cable, the LEDs will light up and start flashing.
- When the device is fully charged, the LEDs stop flashing and the LED indicator lights off.

Caution! Charging time – **approx. 130 minutes.**

**Attention!** Recommended ambient temperature when charging the battery: between 0°C and + 45°C.

## 7.2 Turning on/off

1. To turn on the device, press and hold on the on / off / switch function selector (5) until you hear a beep. The LED will light up with the cleansing / smoothing function (Peeling / lifting). After that you can select the function you prefer by pressing shortly on / off / function switch button.
2. Each subsequent press of the on / off / switch function button changes the functions in sequence: Peeling / lifting, EMS / Pulse, Ion +, Ion-. Every change is signaled by the lighting of the appropriate LED diode.
3. To turn off the device, press and hold the on / off / switch function selector (5) until the illumination is off. Switching off the device will be signaled by a characteristic sound.

## 7.3 Usage

Attention! Before using the device:

Make sure that you do not show any contraindications listed in section **2 Safety instructions**.

Take off your jewelry.

Put a headband on your hair to keep it from falling over your forehead.

Remove make-up / cleanse the skin of the face, neck and décolleté.

Remove the spatula cover.

Disinfect the spatula.

Turn on the device, then select the appropriate function and proceed with the treatment.

### **The device has 4 functions:**

1. Cleansing / smoothing (Peeling / lifting)
2. EMS massage (EMS / Pulse)
3. Positive ionization (Ion +)
4. Negative ionizing (Ion-)


#### **1. Cleansing / smoothing (Peeling / lifting) -**

This function ensures a deep refreshment of the epidermis and the firming of the skin by using ultrasounds and vibrations that remove the stratum corneum and soften the skin.

The best is to be used on damp skin. Wet the skin with warm water or use a facial cleanser.

You can also apply a cleansing gel to create a foam and remove it with the device by putting a spatula close to the skin as shown in the figure below. Then slowly move the spatula against the pores to clean for 3 minutes.

Correct position of the spatula on the skin:

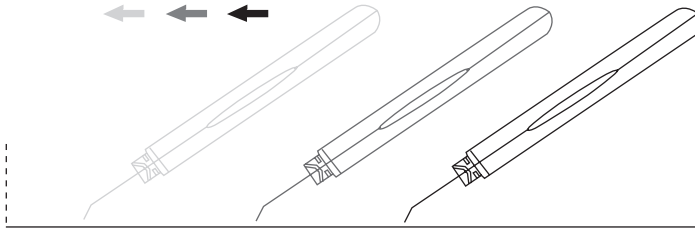Advice:


Start cleansing from the nose. Left side of the nose, bridge, right side of the nose. Always bottom up. Then part of the face from the jaw line to the line under the eye and cheeks, bottom up, on both sides.

The next step is the chin, then the skin over the mouth and forehead, also lead the wheel from the bottom up. After cleansing, wash the skin with clean water or a tonic.

It is good to apply a mask, serum or collagen with elastin after the treatment, but a moisturizing cream is also enough.


2. **EMS / Pulse - EMS massage** -electrical muscle stimulation - restores skin elasticity through the emission of delicate and fast electrical impulses transmitted to the muscles. Electrostimulation treatments on the face strengthen the muscles, restore their natural elasticity, which in turn raises the face oval and guarantees a young and healthy appearance.

Correct position of the spatula:

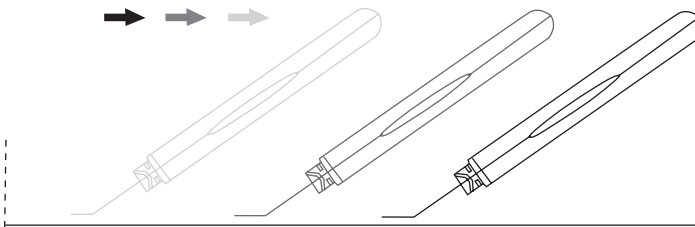
3. **Positive ionization (Ion +)** -deeply cleanses the pores and detoxifies the skin by exporting positive ions. To be used with a toner or cleansing milk.

Correct position of the spatula:

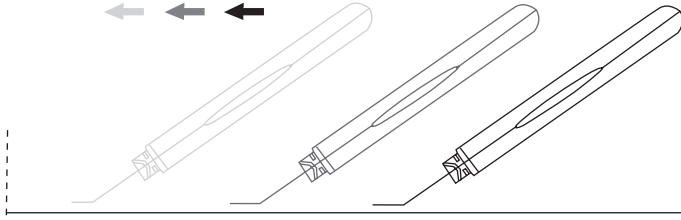4. **Negative ionization (Ion-)** - introducing nutrients deep into the skin.

Correct position of the spatula:

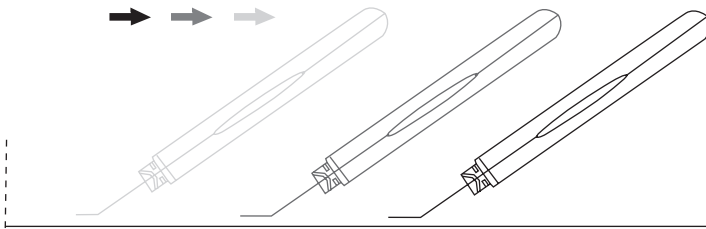Start the treatment by slowly moving the spatula over the skin.  
Do not press the device too hard against the skin.

### Caution!

1. The treatments can be conducted in any order. It is not necessary to maintain the order or use all functions during the treatment.
2. **Caution!** When in use, the angle of the spatula should be operated as shown in the pictures (about 45 degrees).
3. When using, hold the conductive strips on the sides of the device with your hand, so that the ion import and export function can work properly.
4. Don't stay in a fixed place on the skin but move the spatula slowly.
5. If there is inflammation or irritation on the skin, wait until the symptoms disappear before using.
6. Small irritation after using the device is a natural process, but if you notice that it has not disappeared for few hours-do not use the device again and contact your doctor.
7. It is recommended to use the device 2-3 times a week.
8. Do not use it multiple times in one day.
9. Always disinfects the spatula before and after use.

## 8. CLEANING AND MAINTENANCE

Correct and regular cleaning will improve the safety and extend the lifespan of the device.

### **Caution! Turn off the device.**

1. The scraper must be disinfected with any disinfectant or hydrogen peroxide after each use.
2. Leave the parts to dry before using the device again.
3. Do not use any corrosive or abrasive agents for cleaning the device.
4. Do not immerse the device in water or other liquids.

## 9. REPAIR AND MAINTENANCE

The device does not contain any user-serviceable parts. Do not attempt to repair the device on your own. Always have the device repaired by a professional.

## 10. STORAGE AND TRANSPORT

1. Keep the device in a dry and well-ventilated place inaccessible to children.
2. Protect the device from vibration and shock during transport.
3. Store the device with the cover capped on the spatula.

## 11. DISPOSAL

All the packaging materials are 100% recyclable and are labelled as such. The packaging should be disposed of in accordance with local regulations. Keep the packaging materials out of children's reach, as these materials can pose a hazard.

### **Proper disposal of the device:**

1. According to the WEEE Directive (2012/19/EU), the crossed-out whee bin symbol (shown on the side) is used to label all electric and electronic devices requiring segregation.
2. Do not dispose of the spent product with domestic waste: turn it over to an electric and electronic devices collection and recycling point. The crossed-out whee bin symbol placed on the product, instruction manual or package communicates this requirement.
3. The plastics contained in the device can be recycled in accordance with the specific marking. Thanks to recycling, re-use of materials or other forms of utilizing used equipment, you make a significant contribution to the protection of our natural environment.

- Information on electric and electronic devices collection locations is available from local government agencies or from the dealer.
- Used, spent batteries and rechargeable batteries must be discarded in specially marked containers, handed over to special waste collection points or to electrical equipment dealers.

## 12. CE DECLARATION

This device has been designed, manufactured and marketed in compliance with the requirements of the New Approach Directives, therefore, the product has been marked with the CE mark, and a declaration of conformity has been issued for it, which is made available to market regulators.

## 13. WARRANTY

In order to claim a product, you must **deliver the product to the Customer Service Center at any Biedronka store.**

If you have any product-related questions or problems, write us at:

- [infolinia@vershold.com](mailto:infolinia@vershold.com)
- or call us at: +48 667 090 903.

Your opinion is important to us. Please rate our product on our page:

[www.vershold.com/opinie](http://www.vershold.com/opinie)

- The Manufacturer of this product grants a 36-month warranty valid from the purchase date. Should you find any defect, submit a complaint at the place of purchase.
- A defective/damaged product means a product failing to provide the features described in the Instruction Manual due to intrinsic properties of the appliance.
- The Beneficiary under the warranty shall deliver the appliance to the Customer Service Point at a relevant store in order to submit a complaint about the equipment. When asserting your rights under the warranty, you need to present proof of purchase (purchase receipt, invoice) together with a description of the device defect in as much detail as possible, in particular of the external manifestations of the defect.
- The warranty excludes products with physical damage other than that caused by the manufacturer or distributor, including specifically any damage resulting from incorrect usage or force majeure.
- The warranty does not cover products from which the tamper-proof seal has been removed and products to which repairs, alterations or structural modifications have been attempted.

6. The warranty does not cover components subject to natural wear and tear during operation.
7. It is recommended that the complete product be returned in order to facilitate defect verification.
8. This warranty for the sold consumer product does not exclude, restrict or suspend the rights of the buyer related to the regulations concerning implied warranty for defects pursuant to the Civil Code of 23 April 1964 (Journal of Laws Dz.U.2014.121).
9. The warranty shall be applicable within the territory of the Republic of Poland.

**Manufacturer (“Warrantor”):**

VERSHOLD POLAND Sp. z o.o.

ul. Żwirki i Wigury 16 A

02-092 Warszawa, Poland

Made in China