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SPIS TREŚCI

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User Manual
Sonic toothbrush

Table of Contents	
1. INTRODUCTION	14
2. INTENDED USE	14
3. TECHNICAL DATA	14
4. PRINCIPLES OF SAFE USE	15
5. EXPLANATIONS OF THE SYMBOLS	17
6. CONTENTS OF THE SET	18
7. CONSTRUCTION	19
8. USE	19
9. CLEANING AND MAINTENANCE	22
10. TROUBLE SHOOTING	23
11. STORAGE AND TRANSPORT	23
12. WARRANTY AND SERVICE	24
13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS	24
14. MANUFACTURER	24

1. INTRODUCTION

Thank you for choosing our product. We are confident that its use will be a great pleasure for you. Please read the user manual carefully before using the product, in particular the instructions on safe use. Please keep this manual in order to be able to refer to it while using the product.

2. INTENDED USE

The sonic toothbrush is intended exclusively for personal oral cavity care in the household. Do not use the sonic toothbrush commercially, industrially, or in medical or veterinary practice.

3. TECHNICAL DATA

Sonic toothbrush	
Model	AL-010125
Lot number	202109
Number of pulsations	42 000 micro vibrations per min.
Number of programs	4
Program duration	2 minutes with 30 seconds interval
Use time after full charge	ca. 90 days (10 days x 2 minutes x 2 cleaning modes a day)
Ingress protection	IIPX7
Cradle charger power supply	5 V =, 3 W

Built-in battery	Li-ion 3.7 V = 1200 mAh, 4.44 Wh
Accessories	2 interchangeable tips, charging base
Colour versions	White, black

4. PRINCIPLES OF SAFE USE

Read the entire user manual carefully before use.

This appliance may be used by children of at least 8 years of age and people with reduced physical and mental capabilities, or lack of experience and knowledge of the appliance, if supervised; or instructed about safe use so that the hazards involved are understood.

Unattended children must not clean or maintain the unit.

Children must not play with this device.

The device is intended for use only with the included charging base. Do not use other charging bases to charge the device.

The integral power supply cable cannot be replaced. If the power supply is damaged the device becomes unusable.

- Use this device only when it is properly assembled (also applies to the charging base).
- Use the device only as described in these user manual!
- Do not use a damaged device. Do not repair yourself if damaged. Take the device to a professional for repair or deliver it to the place of purchase.
- Do not use this device in locations where it is prohibited to use electrical equipment, or where it may cause harmful interference or other hazards.
- Do not switch on or use the device close to flammable materials.
- Do not use the machine to clean surfaces other than teeth, gums and tongue.
- Keep the device out of the reach of children.

- The manufacturer is not responsible for any problems caused by unauthorized modifications to the equipment.
- When using the device, ensure that there are no objects left near the place where the toothbrush will be placed that could be damaged, broken or knocked over. The power supply base cable should be routed so that it does not come into contact with objects that generate high temperatures. It should also not be subjected to crushing, pulling, or other damage. The charging base cable should not interfere with the movement of persons and pets, as this can lead to tripping, tipping over or injury.
- Avoid exposing the device to very high/low temperatures (below 0°C/32°F and over 40°C/104°F). Plastic parts may deform under extreme temperatures, which may affect the operation of the device or lead to complete damage.
- Do not expose the charging base and its power cord to water or other liquids. This may adversely affect the operation and lifetime of equipment components or cause a short circuit or even electric shock. Do not immerse the brush in water deeper than 1 m or for a period longer than 30 minutes.
- Avoid exposing the battery to very high/low temperatures (below 0°C/32°F and over 45°C/113°F). Extreme temperatures may influence capacity and battery life. Avoid contact with liquids and metal objects as it may cause partial or complete battery damage. Use the battery only for its intended purpose. Do not destroy, impair or throw the battery into flames - it might be dangerous and may cause fire. Worn-out or impaired battery should be placed in a dedicated container. Never open the battery. Please dispose of battery according to the instructions. Disconnect an unused cradle charger from power grid.

CAUTION: Risk of fire and burns. Follow manufacturer's instructions.

5. EXPLANATIONS OF THE SYMBOLS

	<p>The device is marked with a crossed-out rubbish bin, in accordance with the European Directive 2012/19/WE on used electrical and electronic Equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled. Collecting this type of waste in designated locations, and the actual process of their recovery contribute to the protection of natural resources. Proper disposal of used electrical and electronic equipment is beneficial to human health and environment. To obtain information on where and how to dispose of used electrical and electronic equipment in an environmentally friendly manner user should contact their local government office, the waste collection point, or point of sale, where the equipment was purchased.</p>
	<p>In accordance with directive 2013/56/UE on the disposal of batteries and accumulators, this product is marked with the symbol of a crossed-out waste bin. The symbol indicates that the batteries used in this product should not be disposed with regular household waste, but treated according to the legal directives and local regulations. It is not allowed to dispose of batteries and accumulators with unsorted municipal waste. Users of batteries and accumulators must use the available collection points of these items, which allows them to return, recycle and disposal. Within the EU, the collection and recycling of batteries and accumulators is a subject to separate procedures. To learn more about existing procedures of recycling of batteries and accumulators please contact your local office or an institution for waste disposal or landfill.</p>
	<p>The mark placed on the product indicates that the product complies with the basic requirements and other provisions of the so-called "Directive of the 'New Approach' Directive of the European Union.</p>
	<p>The packaging of the device can be recycled in accordance with local regulations.</p>
	<p>The packaging elements should be disposed of in appropriate containers.</p>

		Marking of the container to which the cardboard box and insert should be disposed of – paper
	Symbol for cardboard (packaging material).	
IPX7 (toothbrush)	IPX7 degree of ingress protection, according to which the device is waterproof to a maximum depth of 1 m for no longer than 30 minutes (applies to the lowest part of the device).	

6. CONTENTS OF THE SET

After unpacking, check that all listed parts are in the package and that the device is working properly. Included in the set:

1. Sonic toothbrush
2. Interchangeable tips (2 pcs)
3. Charging base
4. User Manual

If any of the above items is missing or damaged, please contact your dealer.

We advise to keep the container in case of any complaints. Keep the packaging away from children - risk of swallowing and/or suffocation.

You print a version of the user manual with the large font size or download it from the web address below:


<https://mptech.eu/qr/hoffen>

Replaceable brush available for purchases from the manufacturer.

Contact: (+48 71) 71 77 400 or pomoc@myphone.pl

7. CONSTRUCTION

GENERAL APPEARANCE OF THE DEVICE AND ITS ACCESSORIES

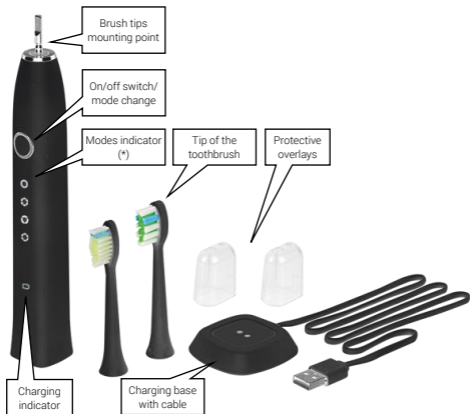*) During operation of the unit, only the current mode indicator is visible.

8. USE

Before and after work, check each component for signs of damage. Damaged components must not be used, as this may cause damage to the user's body.

Do not use the tip for more than 3 months after first use. Teeth or gums might get injured due to excessive use of single brush or bristles wearing off.

MOUNTING AND DISMOUNTING OF BRUSH TIPS

Take out the toothbrush and one of the tips from the packaging. Take off protective cover and put on the tip with the hole on the protruding metal spindle from the brush holder so that the bristles of the tips point in the same direction as the switch and diodes on the brush holder. Press the two parts together in a firm motion. A small space of approx. 0.5 mm will be left between the brush tip and the brush holder. This is the result of the design of the device. Do not press against the brush holder forcing the tip to come into contact with the brush.

To remove the tip, gently pull the tip upwards until it is detached from the brush holder.

CHARGING THE TOOTHBRUSH BATTERY

Keep your hands dry before operating the charging base and charging the battery in the device.

Before charging the battery for the first time, remove the charging base from the package. Unwind the charging base cable gently and place it on a stable horizontal surface. Connect the USB connector of the base power cord to a powered USB adapter (not supplied) with an output of 5V – 0.5A or higher.

Place the bottom of the brush holder in the charging base dock. Make sure that the charging contacts on the underside of the toothbrush holder line up with the contacts on the base. The base has a built-in magnet, which makes the connection to the toothbrush handle more stable. When charging starts, the diode starts to pulse green. In the event of extreme battery discharge, the device may not turn on or might turn off in the middle of a cleaning programme, not completing it, or not turn on at all. When charging is complete, the charging indicator will be continuously illuminated in green.

Do not charge or store the unit below 0°C or above 40°C, otherwise the quality and length of the charge will be affected.

The device may become slightly hot during charging; which is normal.

Try not to discharge the battery completely, as this can lead to permanent damage to the battery and brush.

The built-in battery is not user replaceable. To replace the battery, contact the manufacturer's service department.

TOOTHBRUSH MODES

The toothbrush offers four modes of operation:

- **Cleaning** - strong vibrations
- **Whitening** - medium vibrations
- **Care** - low vibrations
- **Massage** - mixed vibrations

Pressing the switch button again immediately after turning on the device will switch between modes. When all modes have been activated, pressing the switch again will switch the brush off. Operating a mode for approx. 3 seconds activates it and you cannot switch between the modes without switching off the toothbrush or terminating the programme and selecting the programme again.


The toothbrush has a built-in timer and after every 30 seconds it informs you with a short break that the time has passed. The standard washing session provides 4 such intervals, which gives a total of 2 minutes of washing/massage.

TECHNIQUE OF CLEANING WITH A SONIC BRUSH

To use the sonic toothbrush effectively, follow these instructions.


Apply a small dose of your favourite toothpaste to the bristles of the brush. In the case of sonic toothbrushes, the amount of toothpaste may be smaller than in the case of ordinary toothbrushes.

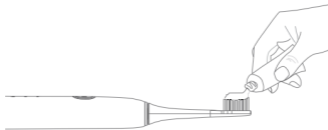The bristles of the toothbrush should be directed at an angle of 45 degrees so that the bristles have the ability to clean the side surfaces of the tooth and gums.

Switch on the device and guide the toothbrush tip slowly, pressing it gently against the tooth surface and touching gums. Single cleaning cycle lasts 4x30s. It is worth dividing the surface of your teeth into four parts so that you can clean each one for 30 seconds. After 2 minutes of operation, the brush will turn off.

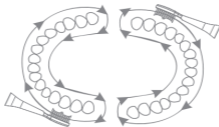During the first dozen or so uses of the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon)

among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. If, after 2 weeks of using the toothbrush in the recommended way gums still bleed, stop using the device and consult a specialist dentist.

9. CLEANING AND MAINTENANCE

CLEANING OF THE EQUIPMENT COMPONENTS

Regular and proper cleaning and maintenance will ensure a long service life and trouble-free operation of the equipment.

Cleaning and maintenance may only be carried out when the appliance is switched off.

It is not allowed to clean any part of the toothbrush in the dishwasher.

Do not immerse the brush in water deeper than 1 m or for a period longer than 30 minutes.

The unit's components can be cleaned under running water at low pressure and temperature. In the case of heavy soiling, mild detergents, such as non-concentrated dishwashing liquids, can be used. Then thoroughly dry the surfaces of the components with a paper or textile towel and allow to dry completely in an airy place. The brushes should be washed under running water and then allowed to dry. From time to time you can wash the tip using a mouthwash.

10. TROUBLE SHOOTING

In case of any problems with the device, please refer to the following tips.

The brush does not switch on	Check that the battery is charged. Perhaps the service life of the battery is over - contact the manufacturer's service department.
Vibrations are too delicate	Check that you are not pressing the brush against the teeth too hard. Check whether the brush tip or brush housing are damaged. Check that the tip is not accidentally too firmly seated on the brush's pin, which may hinder its movement.
Despite the toothbrush has been placed in a charging base, the brush does not charge.	Check that the cable or plug coming out of the charging base is not damaged. Check if the charging base is correctly connected to the power supply. Perhaps the service life of the battery is over - contact the manufacturer's service department.
The gums bleed after using the toothbrush.	Check that you do not press the tip too hard against your gums. Check that you do not press the tip against your gums at a bad angle. <i>During the first dozen or so uses of the gums may bleed slightly (not very frequent phenomenon) or you may feel discomfort (a frequent phenomenon among people who have not used mechanical toothbrushes so far). This usually happens when the gums are not used to the new washing technique and the amount of brush vibrations. If, after 2 weeks of using the toothbrush in the recommended way gums still bleed, stop using the device and consult a specialist dentist.</i>
If the problem persists, contact the manufacturer's service department.	

11. STORAGE AND TRANSPORT

STORAGE OF EQUIPMENT

If you are not going to use the device for a long time, after disconnecting the charging base from the power source, cleaning and drying its components, store it in a dry room away from heat and moisture sources. It is best to place the device in its original packaging, away from children and pets.

TRANSPORT OF THE DEVICE

It is recommended to transport the device in its original packaging. Do not expose the unit to excessive vibration or shock during transportation. During operation, only move the unit when it has been disconnected from the power supply. Do not expose the device to high temperatures or moisture.

12. WARRANTY AND SERVICE

The product (the toothbrush and the charging base) is covered by a 3 year warranty from the date of sale of the product. In case of a defect, the product should be returned to the shop where it was purchased, remembering to provide the proof of purchase. The product should be complete and in its original packaging. Consumables (replaceable brushes) are not covered by the manufacturer's warranty.

13. DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS

mPTech sp. z o.o. hereby declares that the product complies with all relevant provisions of the so-called the 'New Approach' Directive of the European Union.

14. MANUFACTURER

mPTech Sp. z o.o.
Nowogrodzka 31 street
00-511 Warsaw
Poland

Helpline: (+48 71) 71 77 400
E-mail: pomoc@mpotech.eu
Manufactured in China
Lot number: **202109**