HOFFEN

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VERSHOLD POLAND Sp. z o.o.

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
02-092 Warszawa, Polska

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ANALYTICAL SCALE FOR WHOLE-BODY MEASUREMENTS

Model: BFS-2338-B,
BFS-2338-W

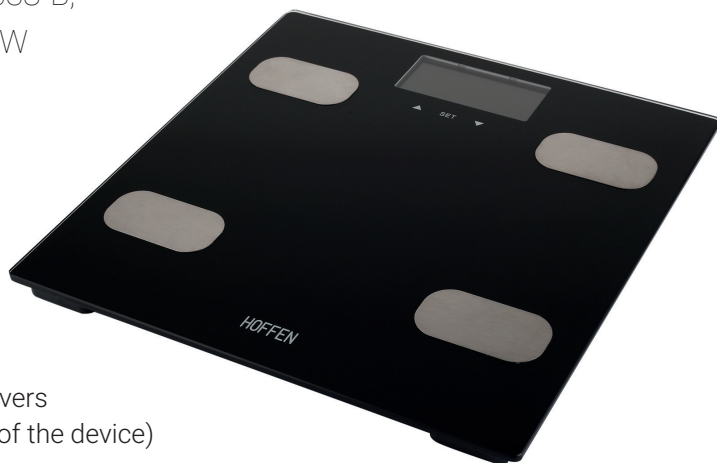(The user's manual covers
various color designs of the device)

TABLE OF CONTENTS

1	INTENDED PURPOSE	22
2	TECHNICAL SPECIFICATIONS	22
3	SAFETY INSTRUCTIONS	22
4	EXPLANATION OF SYMBOLS	24
5	OVERVIEW	25
6	PACKAGE CONTENTS	29
7	USE	29
8	CLEANING AND CARE	35
9	REPAIR	35
10	STORAGE AND TRANSPORTATION	35
11	DISPOSAL	35
12	CE DECLARATION OF CONFORMITY	35
13	WARRANTY	36

1 INTENDED PURPOSE

The analytical scale with whole-body measurements is a device for measuring body weight and percentage measurement of the parameters: body fat, hydration, muscle mass and bone mass. Its simple operation and a special platform design make it a useful household device.

This device is intended for private indoor use only; it may not be used for professional purposes.

The maximum capacity of the scale is 150 kg.

Following the instructions in this manual will ensure safe installation and use of the device.

2 TECHNICAL SPECIFICATIONS

Lot No.	POJM220338
Model	BFS-2338-B, BFS-2338-W
Power supply	3 V  , 2 x battery AAA 1,5 V  (included)
Maximum capacity	150 kg
Memory range	10 users
Height range	from 60 to 240 cm
Age range	recommended for persons between 10 to 85 years old
Weight measurement range	from 2 to 150 kg (if the weight is lower than 5 kg, the display will show C after 1 second, then return to 0,00 kg)
Body fat measurement range	from 3% to 50%
Hydration measurement range	from 25% to 75%
Muscle mass measurement range	from 25% to 75%
Bone mass measurement range	from 1% to 15%

3 SAFETY INSTRUCTIONS

1. The device may be used by individuals with reduced physical, sensory or mental capabilities or by individuals lacking the necessary experience and knowledge, provided that they are supervised or use the device safely, in compliance with the user's manual, and understand the hazards involved.
2. Children must not play with the device.
3. Please read this manual in whole before using the product.
4. **Warning!** The scale is not intended for the following users:

- Children under 10 years of age
 - Persons over age 85
 - Persons with such symptoms as fever, oedema (swelling)
 - Persons with osteoporosis
 - Persons who undergo dialyses regularly
 - Persons with a pacemaker
 - Persons who take cardiovascular medications
 - Pregnant women
 - Persons who train intensively longer than 5 hours a day
 - Persons whose heart rate is lower than 60 bpm.
5. When using electrical devices, especially in the presence of children, basic safety precautions should be adhered to.
 6. The device is intended for private indoor use only.
 7. Do not use the device for any purpose other than its intended use.
 8. Do not use electric devices when taking a bath or shower.
 9. Do not use the device in the vicinity of heat sources (heaters, fire).
 10. Do not place any open flame, for example lit candles, near the device.
 11. Remove the batteries whenever the device is not to be used for any significant length of time or before cleaning.
 12. Do not use any caustic agents for cleaning the device.
 13. Protect the device against penetration by water or other liquids.
 14. Do not immerse the device in water or other liquids.
 15. Protect the device against moisture, dripping or splashing water. Do not place objects filled with water, e.g. vases, on or near the device.
 16. To avoid accidents, do not place the device on a slippery surface.
 17. Do not jump on the scale.
 18. The device is not a toy.
 19. Do not step on the bathroom scale with bare and wet feet – you could slip and fall.
 20. Retain the user's manual and the packaging, if possible.
 21. The scale's maximum capacity is 150 kg.
 22. If you hand this device over to a third party, transfer this user's manual as well.
 23. To obtain an electronic version of the user's manual, visit: <http://instrukcje.vershold.com>

4 EXPLANATION OF SYMBOLS

	<p>Read this manual.</p>
	<p>The product is compliant with the requirements of the applicable EU directives.</p>
	<p>Disposal of waste electrical and electronic equipment – see the DISPOSAL section of this manual.</p>
	<p>Instructions for the disposal of batteries – see the DISPOSAL section. The chemical symbol Pb means that the lead content of the battery exceeds 0.004%.</p>
	<p>Trademark indicating that the manufacturer has made a financial contribution to the development and operation of a packaging materials recovery and recycling system.</p>
	<p>Designation of the packaging material: corrugated cardboard.</p>
	<p>Caution! Slippery surface.</p>
	<p>For indoor use only.</p>
	<p>Segregation sign – Labeling indicating the need for waste segregation.</p>
	<p>Designation of the container to which the waste should go – paper.</p>
	<p>Designation of the relevant waste bin – plastic.</p>



Symbols in monochrome printout only for illustration.

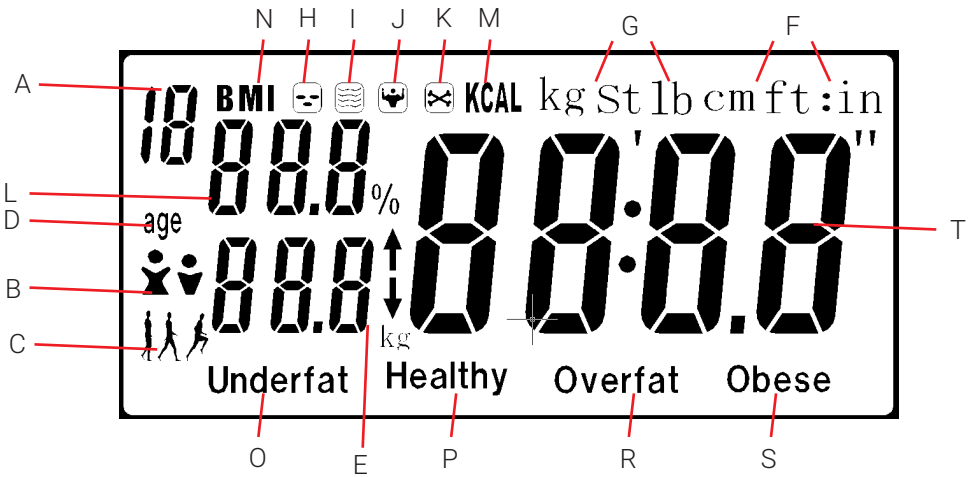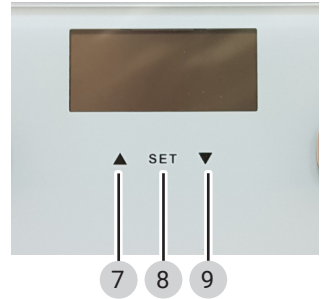5 OVERVIEW

1. Display
2. Top of the scale
3. Metal sensors
4. Battery compartment
5. Battery compartment cover
6. Feet
7. ↑ touch button "up"
8. SET touch button "set"
9. ↓ touch button "down"

	SYMBOL ON THE DISPLAY	MEANING
Selection and setting buttons		
A		User selection (from 1 to 10)
B		Sex selection
		Male
		Female
C		Body composition type
		Sedentary
		Average
		Fit
D	age	Age of user
E		Desired weight – target
F	cm (ft:in)	Height
G	Kg (St:lb)	Unit selection: kilograms, pounds or stones (if you select kilograms, the height unit is provided in centimetres; if you change the weight unit to pounds or stones, height is automatically converted to feet 'and inches')

Measurement display icons

H		Body fat measurement
I		Measurement of water content in the body
J		Muscle mass measurement
K		Bone mass measurement
L	%	Percentage content of fat/muscles/bones/ water in the body
M	KCAL	Daily calorie intake
N	BMI	Body mass index (body weight in KG divided by the square of the height in M)
Description of measurement in relation to the planned weight		
O	Underfat	Underweight
P	Healthy	Healthy – correct weight
R	Overfat	Overweight
S	Obese	Obesity
T		Height/Scale measurement

6 PACKAGE CONTENTS

QTY	Item
1	Bathroom scale with measurement of parameters
2	AAA 1,5 V --- battery
1	User's manual

Open the packaging and carefully take out the product. Check that the set is complete and intact. Check that the plastic parts are not cracked, and the power cord is not damaged. If any part is damaged or missing, contact the dealer and do not use the product. Keep the package or recycle it in accordance with the applicable regulations.

Caution! For the sake of children's safety, do not leave any packaging materials accessible and unattended (plastic bags, cardboard boxes, Styrofoam etc.). Risk of suffocation!

7 USE

Caution! The purchased device may be supplied with batteries installed in the battery compartment and secured with a piece of plastic. Remove the piece of plastic before using the device.

7.1 Installation /replacement

1. Make sure the scale is switched off.
2. To open the battery compartment, lift up its cover.
3. Remove the spent batteries (if present).
4. Insert two new batteries (AAA 1.5 V ---), observing the polarity marking inside the battery compartment.
5. Close the battery compartment cover.

Caution!

- The device is designed to be used with AAA 1.5 V --- batteries.
- The batteries are included.
-  When installing new batteries, observe the markings in the battery compartment.
- Use the type of battery recommended for this device: AAA 1.5 V --- .
- Do not put spent batteries in a rubbish bin, but in special containers for spent batteries.
-  Do not recharge batteries not designed for that purpose (non-rechargeable batteries).

-  Do not recharge batteries not designed for that purpose (non-rechargeable batteries).
-  Never expose batteries to direct heat sources, such as excessive sunlight, heaters or fire.
-   Explosion hazard! Batteries must not be dismantled, put into fire or short-circuited.
- When the device is not in use, always remove the batteries. This will prevent damage to the device in the event that electrolyte leaks from the batteries.
- Remove spent batteries from the device immediately.
- Electrolyte may leak out of spent batteries, causing damage to the device.
- If you touch acid from a battery, rinse your hands under running water. If any acid gets into your eyes, contact a doctor. The acid contained in the battery may cause irritation or burns.
-  Swallowing a battery may have fatal consequences! Keep batteries away from children and pets. Seek medical assistance immediately if a battery is swallowed.

7.2 Weighing

Place the scale on a level, stable and dry surface. In order to ensure a correct reading, do not place the scale on a rug/ carpet. An uneven surface can affect result accuracy.

CAUTION! Do not stand on the scale with wet feet. Always stand on the scale with bare, dry feet (take off your shoes and socks).

Step on the scale with both feet. The scale will switch on automatically. Wait for the measurement. The digits will start flashing before the final result of weighing is displayed.

CAUTION! It is recommended that you stand still during weighing.

The scale will turn off automatically if idle for approximately 30 seconds.

7.3 Changing the weight units

After switching on the scale you can change the weight unit using the “up” button under the display screen. Each press of the button toggles the unit of weight between kg (kilograms), lb (pounds) or st (stones). When you change the weight unit, the unit of measurement changes automatically (if kg, then cm; if pounds or stones, then inches and feet)

7.4 Entering user settings

After switching the scale on, you can enter individual user parameters such as sex, age and height. For this purpose:

- Turn the scale on by touching the “SET” button.
- Press the “SET” button, then select the required parameters by pressing “up” or “down”,

confirming each selection with another press of the “SET” button.

- The parameters are displayed in the following order:
 - A) User 1, 2, 3 and so forth up to a maximum of 10 users.
 - B) Sex
 - C) Body composition type
 - D) Age
 - E) Planned weight - target
 - F) Height
 - G) Unit of measure selection (kg=cm; lb/St=ft:in)
- To access individual user data, press the “SET” button and select the appropriate profile from those flashing on the display.
- Now you may step onto the scale. Remember to keep the scale on an even surface and to step with bare feet on the metal sensors.
- Stay still for a moment. The scale is performing calculations, which is indicated by the subsequently displayed square symbols.
- After the scale completes the calculations, the display will show the data below twice, in the following order:

	Body fat measurement
	Measurement of water content in the body
	Muscle mass measurement
	Bone mass measurement
BMI	Body mass index (body weight in KG divided by the square of the height in M)
KCAL	Daily calorie intake

7.5 Tables of results

Fat content

Age	Female		Male		Result
	Fat content	Water content	Fat content	Water content	
≤30	4.0-20.5	66.0-54.7	4.0-15.5	66.0-58.1	Underweight
	20.6-25.0	54.6-51.6	15.6-20.0	58.0-55.0	Normal weight
	25.1-30.5	51.5-47.8	20.1-24.5	54.9-51.9	Overweight
	30.6-50.0	47.7-34.4	24.6-50.0	51.8-34.4	Obese
>30	4.0-20.5	66.0-54.7	4.0-15.5	66.0-58.1	Underweight
	20.6-25.0	54.6-51.6	15.6-20.0	58.0-55.0	Normal weight
	25.1-30.5	51.5-47.8	20.1-24.5	54.9-51.9	Overweight
	30.6-50.0	47.7-34.4	24.6-50.0	51.8-34.4	Obese

Correct percentage of bone tissue

Bones (kg)			Bones (kg)				
Female	<39	40-60	>60	Male	<54	55-75	>75
	1,7	2,1	2,4		2,4	2,8	3,1

Muscle content

Muscles (%)				Muscles (%)					
Female	Underweight	Normal	Overweight	Obese	Male	Underweight	Normal	Overweight	Obese
	<25	25-27	28-29	>30		<30	31-34	35-38	>39

Body Mass Index

BMI	Underweight	Normal weight	Overweight	Obese
	<18.5	18.5 – 24.9	25 – 29.9	≥30

Recommended daily calorie intake

Table 2. Energy standards for babies, children, and teenagers, determined at the level of the group energy requirement (EER)

Group/age (years)	Body weight (kg)	MJ/daily						kcal/daily					
		Physical activity (PAL)						Physical activity (PAL)					
		low		moderate		high		low		moderate		high	
Babies 0-0.5 0.5-1	6			2,3						550			
	9			3,0						700			
Children 1-3 4-6 7-9	12			4,3	(1,40)					1000	(1,40)		
	19			5,8	(1,50)					1400	(1,50)		
	27	6,3	(1,35)	7,4	(1,60)	8,6	(1,85)	1550	(1,35)	1800	(1,60)	2100	(1,85)
Boys 10-12 13-15 16-18	38	8,7	(1,50)	10,0	(1,75)	11,5	(2,00)	2050	(1,50)	2350	(1,75)	2700	(2,00)
	54	10,9	(1,55)	12,5	(1,80)	14,4	(2,05)	2600	(1,55)	3000	(1,80)	3450	(2,05)
	67	12,4	(1,60)	14,2	(1,85)	16,6	(2,15)	3000	(1,60)	3400	(1,85)	4000	(2,15)
Girls 10-12 13-15 16-18	38	7,6	(1,45)	8,8	(1,70)	10,2	(1,95)	1800	(1,45)	2100	(1,70)	2450	(1,95)
	51	8,7	(1,50)	10,1	(1,75)	11,6	(2,00)	2100	(1,50)	2450	(1,75)	2800	(2,00)
	56	8,9	(1,50)	10,4	(1,75)	11,9	(2,00)	2150	(1,50)	2500	(1,75)	2850	(2,00)

* correct body weight (BMI 18.5-24.9 (kg/m²))

Table 3. Energy standards for males, determined at the level of the group energy requirement (EER)

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Table 4. Energy standards for females, determined at the level of the group energy requirement (EER)

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7.7 Error codes shown on the display

When an error code shows on the display, wait for the device to turn off. Then turn it on again to repeat the operation. If the error code is still displayed, remove the batteries from the battery compartment for 10 seconds.

EEEE	Displayed when the scale's maximum capacity (150 kg) is exceeded.
ERR	The body fat percentage is out of range or the user has not taken off his/her shoes prior to stepping onto the scale.
Lo	Displayed when the batteries are low. Replace the batteries immediately.
8888	Auto-test (when turning on) - serves to verify whether the display is operating properly.
C	When the weight is below 5 kg, after 1 second C will be displayed, then the indication will return to 0.00 kg.

8 CLEANING AND CARE

Remove batteries from the device before you start cleaning it.

Wipe with a slightly damp cloth with some dishwashing agent, then wipe the scale dry.

Do not use any caustic agents for cleaning the device.

Protect the device from penetration by water or other liquids.

Do not immerse the device in water or other liquids.

9 REPAIR

The device does not contain any user-serviceable parts. Do not attempt to repair the device on your own. Always have a specialist perform repairs.

10 STORAGE AND TRANSPORTATION

1. Keep the device in a dry and well-ventilated place inaccessible to children.
2. Protect the device against vibration and shock during transport.
3. Protect the device from falling.

11 DISPOSAL

Keep the packaging materials out of children's reach, as the materials can pose a hazard. The packaging should be disposed of in accordance with local regulations.

Correct disposal of the device:

1. According to the WEEE Directive 2012/19/EU, the crossed-out wheeled bin symbol is used to label all electric and electronic devices requiring segregation.
2. Do not dispose of the spent product with domestic waste: hand it over to an electric and electronic device collection and recycling centre.
3. Hazardous components in electronic and electrical equipment have a negative impact on the environment and human health.
4. Plastics contained in the device can be recycled in accordance with the specific marking. By recycling materials and spent equipment you will help to protect the environment.
5. Information on electric and electronic devices collection centres is available from local government agencies or from the vendor.

6. Worn out and dead batteries must be discarded in specially marked containers, delivered to special waste collection points or to electronic equipment vendors. Lead content of the batteries exceeds 0.004 %.

12 CE DECLARATION OF CONFORMITY

This appliance has been designed, manufactured and marketed in accordance with the requirements of the New Approach Directives; therefore, the product has been marked with the CE mark and a declaration of conformity has been issued for it, which is made available to market regulators.

13 WARRANTY

In order to claim a product, **submit the device to the Customer Service Desk at any Biedronka store.**

Should you have any questions or issues related to product operation or complaint submission, send them to the following e-mail address:

- **infolinia@vershold.com**
- or contact us by phone: +48 667 090 903

Your opinion is important to us. Evaluate our product at: **www.vershold.com/opinie**

1. The Warrantor of this product grants a 36-month warranty valid from the date of purchase. If any defect is found, a complaint about the product should be lodged at the point of purchase.
2. A defective/damaged product means a product failing to provide the features described in the User's Manual due to intrinsic properties of the device.
3. The Beneficiary under the warranty shall deliver the device to the Customer Service Point at a given store in order to submit a complaint about the equipment. When asserting your rights under warranty, you need to present proof of purchase (receipt, invoice) together with a description of the device defect in as much detail as possible, in particular of the external manifestations of the defect.
4. The warranty excludes products with physical damage other than that caused by the manufacturer or distributor, including specifically any damage resulting from incorrect usage or force majeure.
5. The warranty does not cover products from which any tamper-proof seal has been removed, and products to which repairs, alterations or structural modifications have been made or attempted.
6. The warranty does not cover components subject to natural wear and tear during operation.
7. It is recommended that the complete product be submitted together with the complaint so as to facilitate service operations.
8. This warranty for the sold consumer product does not exclude, restrict or suspend any rights of the Buyer stemming from warranty regulations for defects in goods sold pursuant to the Civil Code of 23 April 1964 (Journal of Laws Dz.U.2014.121).
9. Territorial scope of warranty coverage: territory of the Republic of Poland.

Manufacturer (Warrantor):

VERSHOLD POLAND Sp. z o.o.

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02-092 Warsaw, Poland

Made in China

HOFFEN

hoffen.com.pl

#sprytniezhoffen